

CONTENTS PAGE

Client:	Guangzhou Analytical Testing Center_LabTech C/O: Techcomp Limited		
Quotation Number:	19272	Date:	18/07/2025
Job Number:	IN1372	Date:	24/09/2025

- ☒ 1 QUOTATION ✓
- ☒ 2 CUSTOMER ORDER
- ☒ 3 MACHINE SPECIFICATION DOCUMENT_S027_v11
- ☒ 4 LOG SHEET_S028_v3
- ☒ 5 FAULT REPORTING DOCUMENT_S029_v3
- ☒ 6 TEST ACCEPTANCE SPECIFICATION_S030_v9
- ☒ 7 ASSEMBLY AND TEST RECORD_S023_v5
- ☒ 8 SERIAL NUMBER FORM_S053_v4
- ☒ 9 TEST DATA
- ☒ 10 PACKING LIST_S035b_v6
- N/A ☒ 11 INSTALLATION REPORT_S057_v6
- N/A ☒ 12 INSTALLATION ACCEPTANCE AND WARRANTY CERTIFICATE_S071_v3



ADDITIONAL:

Payment Terms:
 50% / 40% / 10%
Other: please specify _____

Penalty Clause Yes / No

Please specify: _____

Blue Book Completion:

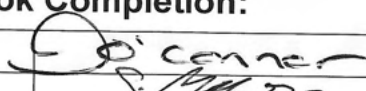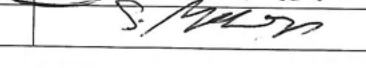Signed		Date	8/4/26
Verified		Date	9/4/26

3 - MACHINE SPECIFICATION

Client:	Guangzhou Analytical Testing Center_LabTech C/O: Techcomp Limited		
Quotation Number:	19272	Date:	18/07/2025
Job Number:	IN1372	Date:	24/09/2025

Ref No.	Product Description	Required
HS2022 (continuous flow)		
20402/1	HS2022 Ref Gas + P Mag	
20402/2	HS2022 Ref Gas + E-Mag	
20402/3	HS2022 Ref Gas + Hd Collector + E-Mag	
20402/4	HS2022 Ref Gas + Hd Collector+E-Mag+Trap Loop	
20501	Hd Collector Assy	***
50101/2	Trapping Loop Assy	
70201-1	70201-1 Analyser spares kit – Included	
EA Units		
10202/11	iso-EArth	
10202/13	iso-EArth+	
10202/14	iso-Earth HT Furnace Upgrade	
10202/15	Iso-Earth (Without Autosampler)	
11802/1	Earth+PC	
11802/2	EArth+Clarity+PC	
70501-1	Consumables for C and N for Combustion Only	
70501-11	Simultaneous CNS Consumables Kit (approx. 10000 analysis)	
70501-3	Sample Prep Kit (Inc.)	
70501-4	GSL Combustion and Gas C & N Start Up Kit	
70501-5	Simultaneous NCS Consumables Kit	
70501-6	C&N Combustion Only Consumables Kit (approx. 5000 samples)	
70501-8	Simultaneous CNS Consumables Kit (approx. 5000 samples)	
70601-1	Sulphur Start Up Kit (Inc. with Sercon EA Systems)	
70601-2	Sulphur Consumables Kit for 1000 Samples	
70601-4	Sulphur Consumables Kit (approx. 5000 Samples)	
70801-2	Pyrolysis Consumables Kit-Low temp for Solids	
71701-1	Sub Microgram Nitrogen Analysis Start Up Kit	
SI2040	PN150 Pneumatic Autosampler with one no.50 pos, carousel	
SI2345	ME Micro Balance	
SI2737/A	Sercon 50 Position Autosampler Assy	
SI2789	Extra 50 Position Carousel	
SI2744	Automatic Liquid Sampler 15 Position	
SI2743	Liquid Sampler Mounting Plate	
Headspace Equilibration		
10202/10	XYZ Gas Sampling Control Box	
10302/1	Aqua System	
11103	µCarbs for CF	
11103-1	µCarbs for DI	
11601	200 Place Dual Heated Block 240 Volt	
11801	Sercon XYZ Gas Autosampler	
71107-1	CF oxygen-18 equilibration kit	
71107-2	DI Oxygen-18 Equilibration kit	
71201-1	CF deuterium equilibration kit	
71201-2	DI Deuterium Equilibration kit	
71801-1	WES-DI start up kit (Inc. with WES-DI)	
GC Units		
10402/2	Orchid module 1500C (Will run Hd 70901-1 C Spares Kit Inc)	
10901	Nitrogen Option (71301-1 N Spares kit Included)	
HP7890/1	Agilent 7890B GC System	
HP7890/2	Agilent 7890B GC with Auto Injector	
HP7890/3	Agilent 7890B GC + Auto Injector + 150 Vial Sample Tray	
SI2561	Nexis GC 2010 Plus	

3 - MACHINE SPECIFICATION

Client:	Guangzhou Analytical Testing Center_LabTech C/O: Techcomp Limited		
Quotation Number:	19272	Date:	18/07/2025
Job Number:	IN1372	Date:	24/09/2025

Ref No.	Product Description	Required
GC Units <i>continued</i>		
SI2564	Nexis GC-2030	
70901-3	GC Pyrolysis H and O Start Up Kit	
70901-4	GC Pyrolysis H Spares Kit	
70901-5	GC Combustion C Spares Kit (approx. 5000 samples)	
70901-6	GC Pyrolysis h Start Up Kit (approx. 5000 samples)	
71301-1	GC N Spares kit	
ABCA2 Units		
10502/1	ABCA2 Sercon XYZ Autosampler	
10602	ABCA2 Gas Box	
20602	Complete Sercon ABCA2 system	
20603	ABCA IR Analyser	
70101-1	ABCA2 spares kit (Inc. With ABCA2)	
70101-2	ABCA2 Spares Kit for AB Analitica Systems	
70101-3	ABCA2 Autosampler Upgrade Kit	
SI1767	CVL505 – Barcode Reader	
SC9973	13C Urea breath test kit for Helicobacter pylori infection	
SC9974	Helicobacter Test kit for children age 3-11 (45mg)	
Integra2 Units		
20802/4	Integra2	
70201-2	Integra2 analyser spares kit - Included	
70401-1	Integra-2 consumables (Inc. With Integra2)	
Trace Gas Systems		
10802/4	CryoPrep	
10802/5	Automated LN refiller	
10802/6	CryoFlex (TOC)	
10802/7	CryoGas	
10802/8	CryoFlex (LA)	
SI2533	Silicon Carbide Vertical Furnace	
SI2731	213nm Laser Ablation System	
SI2733	isoScell 100 sampler chamber	
71001-1	Cryoprep Spares Kit	
71001-5	CryoGas Consumables Kit	
TOC Units		
14001	Thermalox Environmental TOC-IRMS	
14002	Thermalox Environmental TOC-Standalone	
14003	Thermalox Environmental TOC-Standalone No Sampler	
SI2772	CO2 Detector, High sensitivity Option	
SI2773	NOx detector Option	
SI2774	Vacuum pump Option	
SI2775	Stirred Vial rack Option	
72101-1	TOC Start Up Kit	
72101-2	TOC Chem Pack	
72101-3	TOC & TN Instrument with heavy usage - one year	
72101-4	TOC only - economy holding	
GEO 20-22 (dual inlet)		
10702	GEO 20-22 Water Equilibration System (D & 18O)	
21802/1	GEO 20-22	
21802/2	GEO 20-22 + E-Mag	

3 - MACHINE SPECIFICATION

Client:	Guangzhou Analytical Testing Center_LabTech C/O: Techcomp Limited		
Quotation Number:	19272	Date:	18/07/2025
Job Number:	IN1372	Date:	24/09/2025

21802/3	GEO 20-22 + E-Mag +Hd	
21802/4	GEO 20-22 Inlet	
70301-1	GEO start up kit - Included	
70301-3E	Comprehensive GEO 20-22 spares kit for 5 years operation	
HT-EA		
11501/4	HT-EA (With Autosampler)	
11501/5	HT-EA With Liquid Option	
11501/6	HT-EA Stand Alone	
11501/7	HT-EA (Without Autosampler)	
SI2622	CTC PAL3-LSI-L-S	
SI2624	CTC PAL3-Kit-Standalone	
70801-1	Pyrolysis Consumables Kit-High Temp for Solids	
71501-1	HT-EA Pyrolysis Start Up kit O and H	
71501-2	HT-EA Pyrolysis Kit for Water injections	
71501-3	HT-EA Sensitivity Start Up Kit for O & H in small samples	
71501-4	HT-EA Sensitivity Consumables Kit for O & H in small samples	
71501-5	HT-EA Sensitivity Kit for O&H in Small Samples (5000 samples)	
71501-7	HT-EA Kit for Room Temp Equilibration	
71501-8	HT-EA Pyrolysis Kit for Water injections (approx. 1 year kit)	

Ref No.	Product Description	Required
Upgrades		
12101	Evacuation System ***	
12102	2 Ref Gas+ Dilutor Box ***	
20101	Standard Source Assy ***	
20102	HS Source Assy ***	
20301	Standard Collector Assy ***	
30202/1	System Control Upgrade (Pmag) installation not included	Edwards Or Pfeiffer
30202/2	System Controller Upgrade (Emag) installation not included	Edwards Or Pfeiffer
10202/7	Preparation PCB Upgrade	
30202/3	Hijack Box ***	
30302/4	nEXT85 Edwards Pumping System Upgrade Single - installation not included	
30302/5	nEXT85 Edwards Pumping System Upgrade Dual - installation not included	
30302/6	nEXT85 Edwards Pumping System Upgrade Triple- installation not included	
SI2542	ABCA2 Upgrade Kit	
CALLISTOW	Callisto software upgrade ***	
Specials		
11702	LC-I interface	
12202	Leco Interface Unit	
SI2320	Agilent 1260 Series Refractive Index Detection HPLC System	
71401-1	ABCA-G Start Up/Spares Kit	
71401-2	ABCA-G Equilibration Spares Kit - Room temp O only	
71401-3	ABCA-G Equilibration Spares Kit - Room Temp O & H	
71601-1	LC Interface Spares Kit	
ADDITIONS		
PC	Windows 7 based PC System with Callisto Software	
Printer	Printer	
INST1	Installation and start up training at customer's site by Sercon Engineer	
INST2	Installation and start up training at customer's site by Agents Engineer	
WARR	1 year warranty	
WARR2	2 year warranty	
TRAIN	Training in UK for one person for 3 days	
SCONTRACT	Maintenance Service Contract – 1 year	

3 - MACHINE SPECIFICATION

Client:	Guangzhou Analytical Testing Center_LabTech C/O: Techcomp Limited		
Quotation Number:	19272	Date:	18/07/2025
Job Number:	IN1372	Date:	24/09/2025

SPECIFICATION REQUIREMENTS			
Voltage	Volts 50/60Hzs		
Plug Type	SI2510	3 Metre USA/Japan plug to IEC C 13	
	SI2511	2.5 Metre Australian plug to IEC C 13	
	SI2512	3 Metre UK 13 amp Fused plug to IEC C 13	
	SI2523	3metre Schuko (Euro) Plug to IEC C13 Socket	
	SI2514	2.5 metre Chinese mains to IEC C13	
Pallet	Stamped for export heat treated for 30 min @ 56 °C		
	Complies with all Export Requirements		

4 - LOG SHEET


Client:	Guangzhou Analytical Testing Center_LabTech C/O: Techcomp Limited		
Quotation Number:	19272	Date:	18/07/2025
Job Number:	IN1372	Date:	24/09/2025

Date	Action Taken	By Whom
16/2/26	Boxes arrived back	AD
	1 - Caved in top	
	2 - Pallet leg broken	
	3 - Slide off pallet / very damp	
17/2/26	Removed everything from boxes and put MS and Isoearth onto trolleys (1hr)	AD
	Ok removing turbo's, System controller being removed and tested	
	Flight tube visibly bent, and HD headlamp crushed	
	Source and collectors being removed, as well as headlamps.	
	Flight tube being removed from base	
18/2/26	System controller performing ok	AD
19/2/26	All optics and rotary put on Stock 24 for testing	AD
	Emag also put on stock 24 for testing	
20/2/26	Iso earth connected to IN1377 mass spec. Valve connection loose in orange block - reconnected - communicating ok	AD
13/2/26	Source, magnet, collectors and head lamps all responding to tuning on Stock 24	AD
	IsoEarth now has tubes and full helium flow. Pressure test passed.	
	Combustion furnace wouldn't heat, mains wire crimp loose in the strip - Squashed crimp and re-connected.	
	Bracket on furnace base - production will fix	

4 - LOG SHEET

Client:	Guangzhou Analytical Testing Center LabTech C/O: Techcomp Limited		
Quotation Number:	19272	Date:	18/07/2025
Job Number:	IN1372	Date:	24/09/2025

Date	Action Taken	By Whom
23/2/26	Background low 7E ⁻¹¹ 's, with furnaces at temp Banks running overnight	AD
24/2/26	Blanks good. Testing Beet mix Samples. Precision run on - N ₂ failed. Checked ref gas all ok. Adjusted windows, running precision again. Baking out GC overnight	AD
25/2/26	Beet mix precision passes 0.07 CO ₂ and 0.15 N ₂ (manually calculating standard using samples 1, 3, 4, 5, 6) Running another precision overnight	AD
26/2/26	Purged O ₂ and He lines - re run Specs COMS error stopped run. All MS parts on Stock 24 tested and working well, removing to go back on to IN1372 Working on better separation of peaks to include more of N ₂ peaks	AD
27/2/26	N ₂ only passes N ₂ only with dual method - not passing Dual, no CO ₂ source loaded - not passing	AD
2/3/26	Told to stop work on IN1372 units	AD
5/3/26	All parts fitted back to IN1372 Mass spec New / reworked flight tube Fitted pumps, rotary, Sycon, Source, New He collector and Mass Spec on, started pumping down nicely + pumps reach full speed. Emg placed on original workings and left to heat overnight	TI

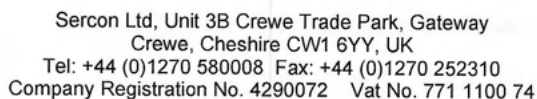

Client:	Guangzhou Analytical Testing Center_LabTech C/O: Techcomp Limited			
Quotation Number:	19272		Date:	18/07/2025
Job Number:	IN1372		Date:	24/09/2025

Date	Action Taken	By Whom
6/3/26	Filament turned on, good emissions Started tuning for CO ₂ + N ₂	Ti
9/3/26	CO ₂ + N ₂ passed ref gas spec's CO ₂ : Prec: 0.03 Lin: 0.19 N ₂ : Prec: 0.04 Lin: 0.17 H ₂ peak detected ISO EARTH connected to main spec. All furnaces heat up as expected Iso crimp for prep unit opened up to ~ 3.0E-6 Pressure test passed with both cold and full temp tubes 1-2 PSI drop over 300 secs N ₂ Blank ~ < 3.0E-10 Dual Beetmix Run passed spec's (PRECISION) CO ₂ : 0.03 N ₂ : 0.13	Ti
10/3/26	ISO EARTH cooled down + turned off HS2022 powered off System passed spec's + ready for packing	Ti
12/3/26	ISO Earth + HS2022 Quality checked Before Packing. System controller rails changed on HS2022. All passed inspection Box 1: HS2022 Box 2: ISO Earth Box 3: Micro Carb, Auto Sampler, Magnet Kit 70201-1, Kit 70501-3, Kit 71107-5, Rotary pump	SW

4 - LOG SHEET

Client:	Guangzhou Analytical Testing Center LabTech C/O: Techcomp Limited		
Quotation Number:	19272	Date:	18/07/2025
Job Number:	IN1372	Date:	24/09/2025

Date	Action Taken	By Whom
5/11/25	Chinese agents setting up mass spec and Isoearth	HC
6/11/25	Source not working great, replaced with low sens one whilst getting repaired	HC
7/11/25	Isoearth Set up and bot, Jumper needed changing on prep board Filament conditioned	HC
10/11/25	tuning ref gas. Pressure test passes	HC
11/11/25	Autosampler on iso needs working as its leaking.	HC
12/11/25	Ref gas for NTC pass.	HC
13/11/25	Blanks look ok, running precision	HC
14/11/25	C passes, N data not good. Checking tubes	HC
17/11/25	H ₂ O trap configuration wrong, trap not being used. Changed to go through trap. Precision passes. Swapped GC column for Sulphur	AD
18/11/25	Swapped source for high Sens one. Conditioned filament	AD
19/11/25	re-tuning ref gas	AD
20/11/25	N ₂ + CO ₂ ref gas pass	
24/11/25	Sulphur peaks found and tuned for Setting up Ucarbs	AD

4 - LOG SHEET

Client:	Guangzhou Analytical Testing Center_LabTech C/O: Techcomp Limited		
Quotation Number:	19272	Date:	18/07/2025
Job Number:	IN1372	Date:	24/09/2025

Date	Action Taken	By Whom
25/11/25	Crimps on ucarbs set. CO-ords for VERSA Set. Full flow through system. needle for flush sampling connected	AD
26/11/25	Pearls coming through on TCD Ucarbs now flushing 5% CO ₂ Flow not coming out of GC completely, tried baking out for a couple hours. GC type not right, swapped for correct one. Flow now fully comes out of GC. Flushing and analysing samples	AD
27/11/25	CO ₂ analysed. C passes, O passes but needs improving so carbonates can pass	AD
28/11/25	CO ₂ ref gas tuned slightly better with good results for dry vials Flushed new vials, left run on over weekend	AD
1/12/25	Iso valve closed on analysis. Re-flushed and ran vials - pass Spas well. Set up acid dose tubing and fittings	Hc/AD
2/12/25	Ran carbonate samples. 500mg carbonate samples meet specification with 0.04% C and 0.08% O. 100mg carbonate samples meet C specification with 0.04% but O fails with 0.21%. 500mg decided as the weight for future analysis	Hc
05/12/25	McArbs system used to run customer samples & calibrate our own calcium carbonate	Hc
08/12/25	HD VEC added to syst con & HD peak shape found	Hc

6 - TEST / ACCEPTANCE SPECIFICATION

Customer:	Guangzhou Analytical Testing Center_LabTech C/O: Techcomp Limited		
Quotation Number:	19272	Date:	18/07/2025
Job Number:	IN1372	Date:	24/09/2025

HS20-22			
Reference Gas Precision (n=5)			
¹³ C	(CO ₂)	0.06	
¹⁵ N		0.06	
¹⁸ O	(CO ₂)	0.06	
³⁴ S	(SO ₂)	0.1	
² H		0.4	

Reference Gas Linearity			
¹³ C	≤0.02‰/nA		
¹⁵ N	≤0.02‰/nA		

IsoEarth & IsoEarth+					
Combustion (n=5)					
¹³ C	100 µg C	0.1		10 µg C	0.8
¹⁵ N	100 µg N	0.15		10 µg N	1.5
³⁴ S	25 µg S	0.3		10 µg S	1.0
Dual	400 µg C	0.1		100 µg N	0.15

Integra2			
¹³ C	100 µg C	0.2	
¹⁵ N	100 µg N	0.3	
Dual	100 µg N & 400 µg C	0.4 & 0.15	
³⁴ S	30 µg S	0.5	
¹⁸ O	100 µg O	0.5	

6 - TEST / ACCEPTANCE SPECIFICATION

Customer:	Guangzhou Analytical Testing Center_LabTech C/O: Techcomp Limited		
Quotation Number:	19272	Date:	18/07/2025
Job Number:	IN1372	Date:	24/09/2025

HT-EA			
Pyrolysis (n=5)			
¹⁸ O	100 µg O or 1 µl water	0.3	
² H	200 µg H or 1 µl water	1.0	

Aqua, µCarbs & GasBox					
Equilibrations (n=5)					
¹⁸ O	0.2 ml water	0.1			
² H	0.2 ml water	2.0			
Carbonates (n=5)					
¹³ C	100 µg CaCO ₃	0.08			
¹⁸ O	100 µg CaCO ₃	0.1			
Gas Analysis (n=5)					
¹³ C	5% CO ₂	0.1		Air CO ₂	0.1
¹⁸ O	5% CO ₂	0.1		Air CO ₂	0.2
² H	0.5 ml H ₂	1.5			
¹⁵ N	0.5 ml N ₂	0.1		100 ppm N ₂ O	0.5
¹³ C	Breath (12 ml)	0.1			
DIC (n=5)					
¹³ C	15ppm (0.5ml)	0.15			

LC-I			
(n=5)			
¹³ C		0.2	

6 - TEST / ACCEPTANCE SPECIFICATION

Customer:	Guangzhou Analytical Testing Center_LabTech C/O: Techcomp Limited		
Quotation Number:	19272	Date:	18/07/2025
Job Number:	IN1372	Date:	24/09/2025

GC Orchid						
Combustion (n=3 for each analyte)						
¹³ C	3 Alkanes	100 ng C	0.2		10 ng C	0.5
¹⁵ N	2 Amines	100 ng N	0.5		14 ng N	1.5
Pyrolysis (n=3 for each analyte)						
² H	3 Alkanes	250 ng H	3.0		25 ng H	10.0
¹⁸ O	3 FAMES	100 ng O	0.5		20 ng O	2.0
Water pyrolysis (n=5)						
² H	0.1 µl water	3.0				
¹⁸ O	0.1 µl water	0.4				

CryoRange						
(n=5)						
¹³ C	Air CO ₂ (12ml)	0.1		Air CH ₄ (125ml)	0.2	
¹⁵ N	Air N ₂ (1ml)	0.1		Air N ₂ O (125ml)	0.4	
¹⁸ O	Air CO ₂ (12ml)	0.15		Air N ₂ O (125ml)	0.7	
¹⁸ O	Air O ₂ (1ml)	0.1				
² H	Air CH ₄	3.0				
¹³ C	TOC (5ppm)	0.4				
¹³ C	Laser					

6 - TEST / ACCEPTANCE SPECIFICATION


Customer:	Guangzhou Analytical Testing Center_LabTech C/O: Techcomp Limited		
Quotation Number:	19272	Date:	18/07/2025
Job Number:	IN1372	Date:	24/09/2025

Geo2			
Internal precision (no change overs)			
² H	200 Bar µl	0.15	
¹³ C	100 Bar µl	0.01	
¹⁵ N	100 Bar µl	0.01	
¹⁸ O	100 Bar µl	0.015	
³⁴ S	100 Bar µl	0.015	
Single cold finger volumes (10 change overs)			
¹³ C	20 Bar µl	0.01	
¹⁸ O	20 Bar µl	0.015	
¹³ C	3 Bar µl	0.02	
¹⁸ O	3 Bar µl	0.03	
CAPS (n=5)			
¹³ C	50 µg CaCO ₃	0.03	
¹⁸ O	50 µg CaCO ₃	0.06	
WES (n=5)			
² H		1.0	
¹⁸ O		0.05	
Manifold (10 tubes)			
² H		0.5	
¹³ C		0.02	
¹⁵ N		0.02	
¹⁸ O		0.03	
³⁴ S		0.03	

6 - TEST / ACCEPTANCE SPECIFICATION

Customer:	Guangzhou Analytical Testing Center_LabTech C/O: Techcomp Limited		
Quotation Number:	19272	Date:	18/07/2025
Job Number:	IN1372	Date:	24/09/2025

Special tests

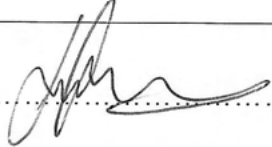Sales Signature: 	Date: 12/12/25
Test Signature: S. Masch	Date: 12/12/25

7: ASSEMBLY AND TEST RECORD

Client:	Guangzhou Analytical Testing Center_LabTech C/O: Techcomp Limited		
Quotation Number:	19272	Date:	18/07/2025
Job Number:	IN1372	Date:	24/09/2025

 PRODUCT: 10202/11-2025/1372

 SERIAL NUMBER: ISO EARTH

 ASSEMBLED:
 (Name) [Signature]

 Date: 27/10/25

 PRE-TEST
 INSPECTION:
 (Production name) [Signature]

 Date: 28/10/25
~~28/10/26~~

 PRE-TEST
 INSPECTION:
 (Test name) [Signature]

 Date: 10/11/25

 SYSTEM TEST:
 (Name) [Signature]

 Date: 12/12/25

 PRE-DESPATCH
 INSPECTION:
 (Name) [Signature]

 Date: 12/12/25

P.A.T. Test results

 In compliance with Low Voltage directive 73/23/EEC
 To Low Voltage standard EN61010 – 1:1993

 Earth Continuity 0.81 Ω

 Insulation >300 $M\Omega$

 Earth Leakage Current 1.00 mA. @ 0.49 kVA.

 Tested by: [Signature]

 Date: 28/10/25

7: ASSEMBLY AND TEST RECORD

Client:	Guangzhou Analytical Testing Center_LabTech C/O: Techcomp Limited		
Quotation Number:	19272	Date:	18/07/2025
Job Number:	IN1372	Date:	24/09/2025


PRODUCT:

MS2022

SERIAL NUMBER:

20402/4-2025/1372

 ASSEMBLED:
 (Name)

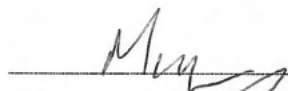

Date:

9/11/25


 PRE-TEST
 INSPECTION:
 (Production name)



Date:

10/11/25

 PRE-TEST
 INSPECTION:
 (Test name)

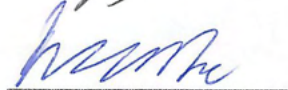

Date:

10/11/25

 SYSTEM TEST:
 (Name)



Date:

12/12/25

 PRE-DESPATCH
 INSPECTION:
 (Name)



Date:

12/12/25

P.A.T. Test results

 In compliance with Low Voltage directive 73/23/EEC
 To Low Voltage standard EN61010 – 1:1993

Earth Continuity

0.056
 Ω

Insulation

>300
 $M\Omega$

Earth Leakage Current


0.30

mA. @

0.24

kVA.

Tested by:

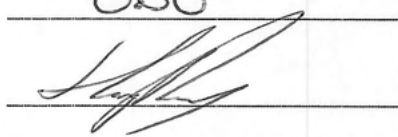

Date:

10/11/25



7: ASSEMBLY AND TEST RECORD

Client:	Guangzhou Analytical Testing Center_LabTech C/O: Techcomp Limited		
Quotation Number:	19272	Date:	18/07/2025
Job Number:	IN1372	Date:	24/09/2025

PRODUCT:

Versa

SERIAL NUMBER:

60101/3-2025/1372

 ASSEMBLED:
 (Name)



Date:

09/10/25

 PRE-TEST
 INSPECTION:
 (Production name)



Date:

7/11/25

 PRE-TEST
 INSPECTION:
 (Test name)



Date:

10/11/25

 SYSTEM TEST:
 (Name)



Date:

12/12/25

 PRE-DESPATCH
 INSPECTION:
 (Name)



Date:

12/12/25

P.A.T. Test results

 In compliance with Low Voltage directive 73/23/EEC
 To Low Voltage standard EN61010 – 1:1993

Earth Continuity

 Ω
NA

Insulation

 $M\Omega$

Earth Leakage Current

mA. @

kVA.

Tested by:

Date:

8 - Serial Numbers Record Document

Client:	Guangzhou Analytical Testing Center_LabTech C/O: Techcomp Limited		
Quotation Number:	19272	Date:	18/07/2025
Job Number:	IN1372	Date:	24/09/2025

Units	S/N
MS/ABCA	20402/4-2025/1372
Integra 2	
Prep Unit 1	10202/11-2025/1372
Prep Unit 2/Misc	11103/2025/1372
Electronics Upgrade	
CETAC	60101/3-2025/1372
Gas Box	
Autosampler	
BCR	
Dual Heater Block	

Pump Kit and Optics	S/N
Turbo pump S	250562847
Turbo pump A	250562848
Rotary Pump	250578697
Voltage	110 / 240
Gauge S	250097503
Gauge A	250452868
Source	20102HS-262
Collector 1	20301-177
Collector 2	20501-082
Magnet	

Parts	S/N
Fan PSU	
24V PSU	

PCBs	S/N
Electronics Box	20250011
Prep PCB	20250011
VFC PCB	25014
Watchdog PCB	
Head Amp PCB 1	25011
Head Amp PCB 2	24611
Electromagnet RP	

Furnaces	S/N
Combustion	
Reduction	
HT-EA	
GC Oven	
Furnace 1	
Furnace 2	

PC	S/N
PC	BKFCF4
Monitor	OMRSGY-TV200-513-4YQB
Printer	

Key

Pumps:	E: Edwards P: Pfeiffer
Furnaces:	V/S: Vecstar E: Elite
Magnet:	E: Electromagnet S: Standard
Turbos & Gauges:	S: Source A: Analyser
Collectors:	Std: Standard S: Sulphur HD: Hydrogen Deuterium

Additional System Info:

Replaced HD headamp - 24012
 Flight tube refurbished - original flight tube used

Packing List

Customer Reference: 42057

Shipper / Exporter:

Sercon Limited
Europa House
Coppicemere Drive
Crewe
CW1 6GZ
UK

Tel: +44 (0)1270 580008
Reference: IN1372

Consignee:

LabTech Corporation, Inc.
No.6 AnQing Street, Area B, Airport Industrial Zone
Shunyi District 101300 Beijing China

Contact Name: Elena Fan
Contact Tel Number: 8610-80486455

BOX: 1

Quantity	Product Code	Product Description	Serial Number: if available
☒ 1	20402/4	HS2022	20402/4-2025/1372
☒ 1		Chinese mains cable	

Packed by:



Date:

15/12/25

Claims for damage or shortages must be made within 3 days

Packing List

Customer Reference: 42057

Shipper / Exporter:

Sercon Limited
Europa House
Coppicemere Drive
Crewe
CW1 6GZ
UK

Tel: +44 (0)1270 580008
Reference: IN1372

Consignee:


LabTech Corporation, Inc.
No.6 AnQing Street, Area B, Airport Industrial Zone
Shunyi District 101300 Beijing China

Contact Name: Elena Fan
Contact Tel Number: 8610-80486455

BOX: 2

Quantity	Product Code	Product Description	Serial Number: if available
☒ 1	10202/11	ISO-EARTH	10202/11-2025/1372
☒ 1		Chinese mains cable	
☒ 1		External cable	

Packed by:

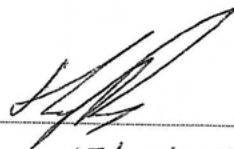

Date:

15/12/25

Claims for damage or shortages must be made within 3 days

Packing List

Customer Reference: 42057

Shipper / Exporter:

Sercon Limited
Europa House
Coppicemere Drive
Crewe
CW1 6GZ
UK

Tel: +44 (0)1270 580008
Reference: IN1372

Consignee:

LabTech Corporation, Inc.
No.6 AnQing Street, Area B, Airport Industrial Zone
Shunyi District 101300 Beijing China

Contact Name: Elena Fan
Contact Tel Number: 8610-80486455

BOX: 3

Quantity	Product Code	Product Description	Serial Number: if available
☒ 1	60101/3	Versa	60101/3-2025/1372
☒ 2		Versa needle	
☒ 1		Versa power supply and Chinese mains cable	
☒ 1		Heater block	

Packed by:



Date:

15/12/25

Claims for damage or shortages must be made within 3 days

Packing List

Customer Reference: 42057

Shipper / Exporter:

Sercon Limited
Europa House
Coppicemere Drive
Crewe
CW1 6GZ
UK

Tel: +44 (0)1270 580008
Reference: IN1372

Consignee:

LabTech Corporation, Inc.
No.6 AnQing Street, Area B, Airport Industrial Zone
Shunyi District 101300 Beijing China

Contact Name: Elena Fan
Contact Tel Number: 8610-80486455

BOX: 4

Quantity	Product Code	Product Description	Serial Number: if available
☒	1	SI2918	PC with Chinese mains cable
☒	1	SI2391	Monitor with Chinese mains cable
☒	1	10202-S	Autosampler
☒	1	SI0685	Electromagnet
☒	1	SI0103	Rotary pump with mains cable
☒	1	11103	uCarbs for CF
☒	1	70201-1	Analyser spares kit
☒	1	70501-3	Sample preparation kit
☒	1	71107-5	uCarbs for CF start up kit
☒	1	70501-1	Start up consumables
☒	1		Software USB (inside PC box)
☒	1		Windows recovery USB (inside PC box)

Packed by:

Date:


15/12/25

Claims for damage or shortages must be made within 3 days