

CONTENTS PAGE

Client:	University of Oxford		
Quotation Number:	19195	Date:	18/06/2025
Job Number:	IN1371	Date:	08/09/2025

- ☒ 1 QUOTATION
- ☒ 2 CUSTOMER ORDER
- ☒ 3 MACHINE SPECIFICATION DOCUMENT_S027_v11
- ☒ 4 LOG SHEET_S028_v3
- ☒ 5 FAULT REPORTING DOCUMENT_S029_v3
- ☒ 6 TEST ACCEPTANCE SPECIFICATION_S030_v9
- ☒ 7 ASSEMBLY AND TEST RECORD_S023_v5
- ☒ 8 SERIAL NUMBER FORM_S053_v4
- ☒ 9 TEST DATA
- ☒ 10 PACKING LIST_S035b_v6
- ☒ 11 INSTALLATION REPORT_S057_v6
- ☒ 12 INSTALLATION ACCEPTANCE AND WARRANTY CERTIFICATE_S071_v3

ADDITIONAL:

Payment Terms:

50% / 50%

Other: please specify

Penalty Clause Yes / No

Please specify:

Blue Book Completion:

Signed	<i>G. O'Connor</i>	Date	26/11/25
Verified	<i>S. Moor</i>	Date	9/4/26

Quotation

To: University of Oxford
 Archaeology Research Laboratory
 Dyson Perrins Building
 South Parks Road
 OXFORD
 OX1 3QY
 VAT No. 195275334

Quotation Date: 18/06/2025
Quotation Reference: 19195
Customer Ref:
Account Ref: UNI009
Quotation Validity: 18/09/2025

Quantity	Product Code	Product Description	Unit Price £	Discount £	Net Price £
Proposal for EAIRMS NERC Funded					
1	20802/4	Integra2 Integrated EA-IRMS, for isotopic and elemental analysis of 13C and 15N, or 34S. Includes standard source, high sensitivity source is available as an optional upgrade. Includes 70201-2 - analyser spares kit - filament, gaskets, pump oil, source tool kit, snoop and ferrules Includes 70501-3 - sample preparation kit - micro spatula, mirror, forceps, sample tray and syringe Includes 70401-1 - Integra2 start up consumables kit Windows based PC	75,000.000	0.00	75,000.00
1	INST	Installation and start up training at Customer's site	7,950.000	0.00	7,950.00
1	WARR2	2 year warranty	0.000	0.00	0.00
1	DEL	Delivery to site	500.000	0.00	500.00

* The price quoted does not include any local taxes or duties
 See separate document for terms and conditions

Payment Details:-

50% with order, 40% on delivery and 10% installation
 Payment is required within 30 DAYS of invoice date by bank transfer.
 All bank transfer charges are the responsibility of the customer.

This quotation is issued in Pound Sterling

Total Net Amount £ 83,450.00

Total VAT Amount £ 16690.00

Order Total £ 100,140.00

Sercon Bank Details:-

National Westminster Bank, 36 High Street,
 Nantwich, Cheshire CW5 5GA U.K.

Sterling: IBAN: GB74NWBK54101736755133
 Euro: IBAN: GB14NWBK60721013806831
 US Dollar: IBAN: GB94NWBK60730113806556
 BIC: NWBKGB2LXXX

121371

STANDARD PURCHASE ORDER: EB2009992

Revision: 0
Document Date: 08-SEP-2025

Sercon Limited
Unit 3B Crewe Trade Park
Gateway
Crewe CW1 6JT
ATTN OF:

Delivery Address
EB04
Research Lab for Archaeology
University of Oxford
Dyson Perrins Building, South Parks Road
Oxford, Oxfordshire OX1 3QY

ATTN OF: Archaeology Institute
ebbilling@arch.ox.ac.uk

Billing Address
School of Archaeology
University of Oxford
1-2 South Parks Road
Oxford, Oxfordshire OX1 3TG

University of Oxford
School of Archaeology
VAT No. GB 125 5067 30

Special Instructions:

EORI No: GB 125 5067 30 000

Supplier Number	Supplier VAT Number	Customer Account No.	Payment Terms	Supplier Contact Phone	Supplier Contact Fax
1024509	GB771110074		30 Days	1270580008	1270252310

Line No.	Part Number/Description	Delivery Date	Qty	UOM	Unit Price GBP	Total GBP	VAT Reference
1	Quotation: 19195 Integrated EA-IRMS	Needed: 22-SEP-2025	1	Each	83,450.00	83,450.00	STANDARD

STANDARD = Standard Rate VAT

Total Excl. VAT:	83,450.00
-------------------------	-----------

Page 1 of 1

Registered Office: Council Sec. University of Oxford Wellington Sq Oxford OX1 2JD United Kingdom

Terms and Conditions

Unless the University has entered into a separate written agreement with you for the supply of goods and services, acceptance of this order implies your acceptance of the University's Standard Conditions of Purchase ([See Terms and Conditions](#)). The University's Standard Conditions of Purchase or, as the case may be, the separate written agreement entered into between the University and you shall govern the agreement for your supply of goods and/or services to the entire exclusion of all other terms and conditions, including, but without limitation, any terms and conditions which you purport to apply under any acknowledgement or confirmation of order, quotation, specification, delivery note, or similar document.

3 - MACHINE SPECIFICATION

Client:	University of Oxford		
Quotation Number:	19195	Date:	18/06/2025
Job Number:	IN1371	Date:	08/09/2025

Ref No.	Product Description	Required
HS2022 (continuous flow)		
20402/1	HS2022 Ref Gas + P Mag	
20402/2	HS2022 Ref Gas + E-Mag	
20402/3	HS2022 Ref Gas + Hd Collector + E-Mag	
20402/4	HS2022 Ref Gas + Hd Collector+E-Mag+Trap Loop	
20501	Hd Collector Assy	***
50101/2	Trapping Loop Assy	
70201-1	70201-1 Analyser spares kit – Included	
EA Units		
10202/11	iso-EArth	
10202/13	iso-EArth+	
10202/14	iso-Earth HT Furnace Upgrade	
10202/15	Iso-Earth (Without Autosampler)	
11802/1	Earth+PC	
11802/2	EArth+Clarity+PC	
70501-1	Consumables for C and N for Combustion Only	
70501-11	Simultaneous CNS Consumables Kit (approx. 10000 analysis)	
70501-3	Sample Prep Kit (Inc.)	
70501-4	GSL Combustion and Gas C & N Start Up Kit	
70501-5	Simultaneous NCS Consumables Kit	
70501-6	C&N Combustion Only Consumables Kit (approx. 5000 samples)	
70501-8	Simultaneous CNS Consumables Kit (approx. 5000 samples)	
70601-1	Sulphur Start Up Kit (Inc. with Sercon EA Systems)	
70601-2	Sulphur Consumables Kit for 1000 Samples	
70601-4	Sulphur Consumables Kit (approx. 5000 Samples)	
70801-2	Pyrolysis Consumables Kit-Low temp for Solids	
71701-1	Sub Microgram Nitrogen Analysis Start Up Kit	
SI2040	PN150 Pneumatic Autosampler with one no.50 pos, carousel	
SI2345	ME Micro Balance	
SI2737/A	Sercon 50 Position Autosampler Assy	
SI2789	Extra 50 Position Carousel	
SI2744	Automatic Liquid Sampler 15 Position	
SI2743	Liquid Sampler Mounting Plate	
Headspace Equilibration		
10202/10	XYZ Gas Sampling Control Box	
10302/1	Aqua System	
11103	µCarbs for CF	
11103-1	µCarbs for DI	
11601	200 Place Dual Heated Block 240 Volt	
11801	Sercon XYZ Gas Autosampler	
71107-1	CF oxygen-18 equilibration kit	
71107-2	DI Oxygen-18 Equilibration kit	
71201-1	CF deuterium equilibration kit	
71201-2	DI Deuterium Equilibration kit	
71801-1	WES-DI start up kit (Inc. with WES-DI)	
GC Units		
10402/2	Orchid module 1500C (Will run Hd 70901-1 C Spares Kit Inc)	
10901	Nitrogen Option (71301-1 N Spares kit Included)	
HP7890/1	Agilent 7890B GC System	
HP7890/2	Agilent 7890B GC with Auto Injector	
HP7890/3	Agilent 7890B GC + Auto Injector + 150 Vial Sample Tray	
SI2561	Nexis GC 2010 Plus	

Ref No.	Product Description	Required
---------	---------------------	----------

S027_v11

Page 1 of 6

23/10/2017

*** When ordered alone IN no. / Blue book NOT required

3 - MACHINE SPECIFICATION

Client:	University of Oxford		
Quotation Number:	19195	Date:	18/06/2025
Job Number:	IN1371	Date:	08/09/2025

GC Units <i>continued</i>			
SI2564	Nexis GC-2030		
70901-3	GC Pyrolysis H and O Start Up Kit		
70901-4	GC Pyrolysis H Spares Kit		
70901-5	GC Combustion C Spares Kit (approx. 5000 samples)		
70901-6	GC Pyrolysis h Start Up Kit (approx. 5000 samples)		
71301-1	GC N Spares kit		
ABCA2 Units			
10502/1	ABCA2 Sercon XYZ Autosampler		
10602	ABCA2 Gas Box		
20602	Complete Sercon ABCA2 system		
20603	ABCA IR Analyser		
70101-1	ABCA2 spares kit (<i>Inc. With ABCA2</i>)		
70101-2	ABCA2 Spares Kit for AB Analitica Systems		
70101-3	ABCA2 Autosampler Upgrade Kit		
SI1767	CVL505 – Barcode Reader		
SC9973	13C Urea breath test kit for Helicobacter pylori infection		
SC9974	Helicobacter Test kit for children age 3-11 (45mg)		
Integra2 Units			
20802/4	Integra2		
70201-2	Integra2 analyser spares kit - Included		1
70401-1	Integra-2 consumables (<i>Inc. With Integra2</i>)		1
Trace Gas Systems			
10802/4	CryoPrep		
10802/5	Automated LN refiller		
10802/6	CryoFlex (TOC)		
10802/7	CryoGas		
10802/8	CryoFlex (LA)		
SI2533	Silicon Carbide Vertical Furnace		
SI2731	213nm Laser Ablation System		
SI2733	isoScell 100 sampler chamber		
71001-1	Cryoprep Spares Kit		
71001-5	CryoGas Consumables Kit		
TOC Units			
14001	Thermalox Environmental TOC-IRMS		
14002	Thermalox Environmental TOC-Standalone		
14003	Thermalox Environmental TOC-Standalone No Sampler		
SI2772	CO2 Detector, High sensitivity Option		
SI2773	NOx detector Option		
SI2774	Vacuum pump Option		
SI2775	Stirred Vial rack Option		
72101-1	TOC Start Up Kit		
72101-2	TOC Chem Pack		
72101-3	TOC & TN Instrument with heavy usage - one year		
72101-4	TOC only - economy holding		
GEO 20-22 (dual inlet)			
10702	GEO 20-22 Water Equilibration System (D & 18O)		
21802/1	GEO 20-22		
21802/2	GEO 20-22 + E-Mag		
21802/3	GEO 20-22 + E-Mag +Hd		
21802/4	GEO 20-22 Inlet		

3 - MACHINE SPECIFICATION

Client:	University of Oxford		
Quotation Number:	19195	Date:	18/06/2025
Job Number:	IN1371	Date:	08/09/2025

70301-1	GEO start up kit - Included	
70301-3E	Comprehensive GEO 20-22 spares kit for 5 years operation	
HT-EA		
11501/4	HT-EA (With Autosampler)	
11501/5	HT-EA With Liquid Option	
11501/6	HT-EA Stand Alone	
11501/7	HT-EA (Without Autosampler)	
SI2622	CTC PAL3-LSI-L-S	
SI2624	CTC PAL3-Kit-Standalone	
70801-1	Pyrolysis Consumables Kit-High Temp for Solids	
71501-1	HT-EA Pyrolysis Start Up kit O and H	
71501-2	HT-EA Pyrolysis Kit for Water injections	
71501-3	HT-EA Sensitivity Start Up Kit for O & H in small samples	
71501-4	HT-EA Sensitivity Consumables Kit for O & H in small samples	
71501-5	HT-EA Sensitivity Kit for O&H in Small Samples (5000 samples)	
71501-7	HT-EA Kit for Room Temp Equilibration	
71501-8	HT-EA Pyrolysis Kit for Water injections (approx. 1 year kit)	

Ref No.	Product Description	Required
Upgrades		
12101	Evacuation System ***	
12102	2 Ref Gas+ Dilutor Box ***	
20101	Standard Source Assy ***	
20102	HS Source Assy ***	
20301	Standard Collector Assy ***	
30202/1	System Control Upgrade (Pmag) installation not included	Edwards Or Pfeiffer
30202/2	System Controller Upgrade (Emag) installation not included	Edwards Or Pfeiffer
10202/7	Preparation PCB Upgrade	
30202/3	Hijack Box ***	
30302/4	nEXT85 Edwards Pumping System Upgrade Single - installation not included	
30302/5	nEXT85 Edwards Pumping System Upgrade Dual - installation not included	
30302/6	nEXT85 Edwards Pumping System Upgrade Triple- installation not included	
SI2542	ABCA2 Upgrade Kit	
CALLISTOW	Callisto software upgrade ***	
Specials		
11702	LC-I interface	
12202	Leco Interface Unit	
SI2320	Agilent 1260 Series Refractive Index Detection HPLC System	
71401-1	ABCA-G Start Up/Spares Kit	
71401-2	ABCA-G Equilibration Spares Kit - Room temp O only	
71401-3	ABCA-G Equilibration Spares Kit - Room Temp O & H	
71601-1	LC Interface Spares Kit	
ADDITIONS		
PC	Windows 7 based PC System with Callisto Software	
Printer	Printer	
INST1	Installation and start up training at customer's site by Sercon Engineer	
INST2	Installation and start up training at customer's site by Agents Engineer	
WARR	1 year warranty	
WARR2	2 year warranty	
TRAIN	Training in UK for one person for 3 days	
SCONTRACT	Maintenance Service Contract – 1 year	
SPECIFICATION REQUIREMENTS		
Voltage	Volts 50/60Hz	

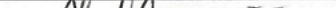
3 - MACHINE SPECIFICATION

Client:	University of Oxford		
Quotation Number:	19195	Date:	18/06/2025
Job Number:	IN1371	Date:	08/09/2025

Plug Type	SI2510	3 Metre USA/Japan plug to IEC C 13	
	SI2511	2.5 Metre Australian plug to IEC C 13	
	SI2512	3 Metre UK 13 amp Fused plug to IEC C 13	
	SI2523	3metre Schuko (Euro) Plug to IEC C13 Socket	4
	SI2514	2.5 metre Chinese mains to IEC C13	
Pallet	Stamped for export heat treated for 30 min @ 56 °C		
	Complies with all Export Requirements		

3 - MACHINE SPECIFICATION

Client:	University of Oxford		
Quotation Number:	19195	Date:	18/06/2025
Job Number:	IN1371	Date:	08/09/2025

Approved by Sales: (Name)		Date: 20/10/25
Approved by Tech.: (Name)		Date: 20/10/25

4 - LOG SHEET

Client:	University of Oxford		
Quotation Number:	19195	Date:	18/06/2025
Job Number:	IN1371	Date:	08/09/2025

Date	Action Taken	By Whom
23/5	Flight tube not sitting flush on top plate. 3-5 mm gap under turbo flange to top.	CS
25/5	Flight tube changed. Holes not aligned, top drilled to fit flight tube hole locations.	CS
5/6	PUMPED DOWN WITH SOURCE + 6 THOU NO COLLECTION FITTED.	PH.
9/7	PC + software set up Gases and connections done Pump on Software v.50.0.3 installed due to EPL issue Crimps set up Sampler leak tight	TI
10/7	Filament conditioned Tubes made and fitted Investigated the gauge issue - all good readbacks okay	TI
11/7	Leak checked system flow ~ 58 ml/min Overnight vacuum 4.2×10^{-8} Oven + GC on - heating up Isolation valve crimp vacuum ~ 4.0×10^{-6}	TI
14/7	Source flashing over, replaced for another - looks ok	HC
15/7	No beams detected	
16/7	Swapped System controller - no beams, Put working head amp on - no beams Put working collector on - no beams	AD/HC

Date	Action Taken	By Whom
------	--------------	---------

4 - LOG SHEET

Client:	University of Oxford		
Quotation Number:	19195	Date:	18/06/2025
Job Number:	101371	Date:	08/09/2025

17/7	Took source out, gas feedthrough knocked out. Put back in, pumped down. CO ₂ peak shape found but not great, flushing over in N ₂ . Put original source on test rig to see how it behaves	AD
18 Jul 25	System flash over. Removed source. filament damaged + swart on IR plate. Fitted new filament. Also adjusted focus magnet position as too close to trap	DJH
18 July 25	System pumping down over weekend	DJH
11.8.25	Unit found pumping down after source mals. Satisfactory tuning for both N ₂ + CO ₂ + 2/1 stability No flash overs Running 400/100 best mix precision overnight.	gpk
12.8.	Good specs for both gas species. Att. Unable to re-process the data file in Re-Processor Software isn't saving the file in Re-Processor.	RK
1-2/19	Fresh tubes high N ₂ background external water trap used + GC baked overnight Connected system to He tank	Ti
Date	Action Taken	By Whom

5 - FAULT REPORTING

Client:	University of Oxford		
Quotation Number:	19195	Date:	18/06/2025
Job Number:	IN1371	Date:	08/09/2025

Fault(s)	Action Taken	By Whom	Total Time Taken to Resolve
18-09-25	Compressed air tubes wrong way round for Autosampler. tubes changed internally	HC	5 minutes
18-09-25	4 Sources & 4 Bakeouts needed for a stable N_2 run without glushovers. Either sources were consistently being made poorly, or flight tube was so unclean it was damaging sources.	HC	5 Days
18-09-25	3 Autosamplers & a replacement of central PTFE slider (back to original white block) needed to stop leaking whilst still turning. New brown sliders found to be too hard & don't compress to stop the leak.	HC	3 days
18-09-25	Leak on trapping loop & O_2 pressure regulator causing high blanks / background	HC	2 hours

6 - TEST / ACCEPTANCE SPECIFICATION

Customer:	University of Oxford		
Quotation Number:	19195	Date:	18/06/2025
Job Number:	ID1371	Date:	08/09/2025

HS20-22			
Reference Gas Precision (n=5) x5 Ref gas pulses of 5 seconds to be used.			
¹³ C	(CO ₂)	0.06	
¹⁵ N		0.06	
¹⁸ O	(CO ₂)	0.06	
³⁴ S	(SO ₂)	0.1	
² H		0.4	

Reference Gas Linearity ref gas pulses ranging from 1-second trough to 10 seconds			
¹³ C	≤0.02‰/nA STD Deviation <0.30		
¹⁵ N	≤0.02‰/nA STD Deviation <0.30		

IsoEarth & PyroEarth					
Combustion (n=5) (For NCS use 0.3mg of Cystine) (Pyrolysis use silver caps) 0.4mg of Cystine for ¹⁸ O and 4 mg of Cystine for H ₂					
¹³ C	100 µg C	0.1		10 µg C	0.8
¹⁵ N	100 µg N	0.15		10 µg N	1.5
³⁴ S	25 µg S	0.3		10 µg S	1.0
Dual	400 µg C	0.1		100 µg N	0.15
¹⁸ O	100 µg or 1 µl water	0.3			
² H	200 µg H or 1 µl water	1.0			

6 - TEST / ACCEPTANCE SPECIFICATION

Customer:	University of Oxford		
Quotation Number:	19195	Date:	18/06/2025
Job Number:	IN1371	Date:	08/09/2025

EArth			
Solid specs to be completed using 2mg methionine calibration against 1,2,3mg (n=3) cystine. Combustion (n=5)			
		Accuracy (avg)	Precision (std dev)
¹³ C		+/- 0.4%	0.2
¹⁵ N		+/- 0.4%	0.2
² H		+/- 0.4%	0.2
³⁴ S		+/- 0.4%	0.2
¹⁸ O	Pyrolysis	+/- 0.4%	0.2
Liquid injections n=5 using 5µl methanol calibrated against 2,4,6µl injections of ethanol			
		Accuracy (avg)	Precision (std dev)
¹⁸ O	Pyrolysis	+/- 0.4%	0.2

Integra2			
Specs to be completed using AS beet mix 400C/100N for precision & 250C/100N for linearity			
¹³ C	100 µg C	0.2	
¹⁵ N	100 µg N	0.3	
Dual	100 µg N & 400 µg C	0.4 & 0.15	0.03% N 0.01% C
³⁴ S	30 µg S	0.5	
¹⁸ O	100 µg O	0.5	
Dual	100 µg N & 250 µg C	0.5	0.12% N 0.30% C
HT-EA (Pyrolysis use silver caps) 0.4mg of Cystine for ¹⁸ O and 4 mg of Cystine for H ₂			
Pyrolysis (n=5)			
¹⁸ O	100 µg O or 1 µl water	0.3	
² H	200 µg H or 1 µl water	1.0	

6 - TEST / ACCEPTANCE SPECIFICATION

Customer:	University of Oxford		
Quotation Number:	19195	Date:	18/06/2025
Job Number:	IN1371	Date:	08/09/2025

Aqua, μ Carbs & GasBox					
Equilibrations (n=5)					
^{18}O	0.2 ml water	0.1			
^2H	0.2 ml water	2.0			
Carbonates (n=5)					
^{13}C	100 μg CaCO_3	0.08			
^{18}O	100 μg CaCO_3	0.1			
Gas Analysis (n=5)					
^{13}C	5% CO_2	0.1		Air CO_2	0.1
^{18}O	5% CO_2	0.1		Air CO_2	0.2
^2H	0.5 ml H_2	1.5			
^{15}N	0.5 ml N_2	0.1		100 ppm N_2O	0.5
^{13}C	Breath (12 ml)	0.1			
DIC (n=5)					
^{13}C	15ppm (0.5ml)	0.15			

LC-I					
(n=5)					
^{13}C		0.2			

GC Orchid					
Combustion (n=3 for each analyte)					
^{13}C	3 Alkanes	100 ng C	0.2	10 ng C	0.5
^{15}N	2 Amines	100 ng N	0.5	14 ng N	1.5

6 - TEST / ACCEPTANCE SPECIFICATION

Customer:	University of Oxford		
Quotation Number:	19195	Date:	18/06/2025
Job Number:	ID1371	Date:	08/09/2025

Pyrolysis (n=3 for each analyte)						
² H	3 Alkanes	250 ng H	3.0		25 ng H	10.0
¹⁸ O	3 FAMES	100 ng O	0.5		20 ng O	2.0
Water pyrolysis (n=5)						
² H	0.1 µl water	3.0				
¹⁸ O	0.1 µl water	0.4				

CryoRange						
(n=5)						
¹³ C	Air CO ₂ (12ml)	0.1		Air CH ₄ (125ml)	0.2	
¹⁵ N	Air N ₂ (1ml)	0.1		Air N ₂ O (125ml)	0.4	
¹⁸ O	Air CO ₂ (12ml)	0.15		Air N ₂ O (125ml)	0.7	
¹⁸ O	Air O ₂ (1ml)	0.1				
² H	Air CH ₄	3.0				
¹³ C	TOC (5ppm)	0.4				
¹³ C	Laser					

Geo2						
Internal precision (no change overs)						
² H	200 Bar µl	0.15				
¹³ C	100 Bar µl	0.01				
¹⁵ N	100 Bar µl	0.01				
¹⁸ O	100 Bar µl	0.015				
³⁴ S	100 Bar µl	0.015				
Single cold finger volumes (10 change overs)						

6 - TEST / ACCEPTANCE SPECIFICATION

Customer:	University of Oxford		
Quotation Number:	19195	Date:	18/06/2025
Job Number:	IN1371	Date:	08/09/2025

¹³ C	20 Bar µl	0.01	
¹⁸ O	20 Bar µl	0.015	
¹³ C	3 Bar µl	0.02	
¹⁸ O	3 Bar µl	0.03	
CAPS (n=5)			
¹³ C	50 µg CaCO ₃	0.03	
¹⁸ O	50 µg CaCO ₃	0.06	
WES (n=5)			
² H		1.0	
¹⁸ O		0.05	
Manifold (10 tubes)			
² H		0.5	
¹³ C		0.02	
¹⁵ N		0.02	
¹⁸ O		0.03	
³⁴ S		0.03	

Special tests			

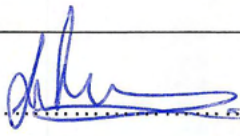
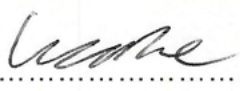


Sercon Ltd, Europa House, Coppicemere Drive,
Crewe, Cheshire CW1 6GZ, UK
Tel: +44 (0)1270 580008
Company Registration No. 4290072 Vat No. 771 1100 74

Test acceptance - 6

6 - TEST / ACCEPTANCE SPECIFICATION

Customer:	University of Oxford		
Quotation Number:	19195	Date:	18/06/2025
Job Number:	IN1371	Date:	08/09/2025

Sales Signature: 	Date: 08/09/25
Test Signature: 	Date: 10/10/25

7: ASSEMBLY AND TEST RECORD

Client:	University of Oxford		
Quotation Number:	19195	Date:	18/08/2025
Job Number:	IN1371	Date:	08/09/2025

PRODUCT:

INTEGRA

SERIAL NUMBER:

20802/4 2025/1371

ASSEMBLED:
(Name)

[Signature]

Date:

4/6/25

PRE-TEST
INSPECTION:
(Production name)

P. Hayes

Date:

11/6/25

PRE-TEST
INSPECTION:
(Test name)

Hayes

Date:

11/06/25

SYSTEM TEST:
(Name)

Hayes

Date:

19/09/25

PRE-DESPATCH
INSPECTION:
(Name)

[Signature]

Date:

23/09/25

P.A.T. Test results

In compliance with Low Voltage directive 73/23/EEC
 To Low Voltage standard EN61010 – 1:1993

Earth Continuity

0.081

Ω

Insulation

7300

MΩ

Earth Leakage Current

0.90

mA. @

40.05

kVA.

Tested by:

[Signature]

Date:

4/6/25

8 - Serial Numbers Record Document

Client:	University of Oxford		
Quotation Number:	19195	Date:	18/06/2025
Job Number:	IN1371	Date:	08/09/2025

Units	S/N
MS/ABCA	20802/4-2025/1371
Integra 2	
Prep Unit 1	
Prep Unit 2/Misc	
Electronics Upgrade	
CETAC	
Gas Box	
Autosampler	42714
BCR	
Dual Heater Block	

Pump Kit and Optics	S/N
Turbo pump S	
Turbo pump A	
Rotary Pump	250175026
Voltage	110 / 240
Gauge S	250067517
Gauge A	
Source	20101HS-257
Collector 1	25010
Collector 2	
Magnet	18-7-28 ⑧

Parts	S/N
Fan PSU	
24V PSU	

PCBs	S/N
Electronics Box	20250005
Prep PCB	
VFC PCB	
Watchdog PCB	
Head Amp PCB 1	
Head Amp PCB 2	
Electromagnet RP	

Furnaces	S/N
Combustion	
Reduction	
HT-EA	
GC Oven	
Furnace 1	
Furnace 2	

PC	S/N
PC	F120154
Monitor	HBWJJ74
Printer	

Key

Pumps:	E: Edwards P: Pfeiffer
Furnaces:	V/S: Vecstar E: Elite
Magnet:	E: Electromagnet S: Standard
Turbos & Gauges:	S: Source A: Analyser
Collectors:	Std: Standard S: Sulphur HD: Hydrogen Deuterium

Additional System Info: