

CONTENTS PAGE

Client:	Central South Forestry University C/O: LabTech, China		
Quotation Number:	18754	Date:	07/03/2025
Job Number:	IN1362	Date:	16/06/2025

- ☒ 1 QUOTATION
- ☒ 2 CUSTOMER ORDER
- ☒ 3 MACHINE SPECIFICATION DOCUMENT_S027_v11
- ☒ 4 LOG SHEET_S028_v3
- ☒ 5 FAULT REPORTING DOCUMENT_S029_v3
- ☒ 6 TEST ACCEPTANCE SPECIFICATION_S030_v9
- ☒ 7 ASSEMBLY AND TEST RECORD_S023_v5
- ☒ 8 SERIAL NUMBER FORM_S053_v4
- ☒ 9 TEST DATA
- ☒ 10 PACKING LIST_S035b_v6
- ☒ 11 INSTALLATION REPORT_S057_v6 *plus Agent install.*
- ☒ 12 INSTALLATION ACCEPTANCE AND WARRANTY CERTIFICATE_S071_v3 *plus Doc. 21/9/25*

ADDITIONAL:

Payment Terms:

50% / 50%

Other: please specify

Penalty Clause Yes / No

Please specify:

Blue Book Completion:

Signed	<i>[Signature]</i>	Date	02/09/25
Verified	<i>S. Melt</i>	Date	03/9/25

Quotation

To: LabTech Corporation, Inc
 No.6 Anqing Street, Area B
 Airport Industrial Park
 Shunyi District
 Beijing, 101300
 P. R. China
 VAT No. 91 110 113 799 019 7

Quotation Date: 07/03/2025
Quotation Reference: 18754
Customer Ref:
Account Ref: LAB014
Quotation Validity: 30/06/2025

Quantity	Product Code	Product Description	Unit Price £	Discount £	Net Price £
Agent Pricing					
1	20802/4	Integra2 Integrated EA-IRMS, for isotopic and elemental analysis of 13C and 15N, or 34S. Includes standard source, high sensitivity source is available as an optional upgrade. Includes 70201-2 - analyser spares kit - filament, gaskets, pump oil, source tool kit, snoop and ferrules Includes 70501-3 - sample preparation kit - micro spatula, mirror, forceps, sample tray and syringe Includes 70401-1 - Integra2 start up consumables kit Windows based PC	65,000.000	0.00	65,000.00
1	70601-3	Sulphur option for Integra Dangerous goods in excepted quantities	3,000.000	0.00	3,000.00
1	70801-4	Oxygen Upgrade Kit for Integra2- Low Temp	880.000	0.00	880.00
1	CON001	Customer spec pc Intel Core i5, Windows 10 ,32GB RAM,512GB SSD 23"Monitor	500.000	0.00	500.00
1	WARR2	2 year warranty	0.000	0.00	0.00

* The price quoted does not include any local taxes or duties
 See separate document for terms and conditions

Payment Details:-

50% with order, 40% on delivery and 10% installation
 Payment is required within 30 DAYS of invoice date by bank transfer.
 All bank transfer charges are the responsibility of the customer.

This quotation is issued in Pound Sterling

Total Net Amount £ 69,202.00
Total VAT Amount £ 0.00
Order Total £ 69,202.00

Sercon Bank Details:-

National Westminster Bank, 36 High Street,
 Nantwich, Cheshire CW5 5GA U.K.

Sterling: IBAN: GB74NWBK54101736755133
 Euro: IBAN: GB14NWBK60721013806831
 US Dollar: IBAN: GB94NWBK60730113806556
 BIC: NWBKGB2LXXX

PURCHASE ORDER

Purchase Order: 20250500013

Date: June 16, 2025

Page: 1

To: Sercon Limited

Unit 3B Crewe Trade Park, Gateway, Crewe, CW1 6JT U.K.

Tel: +44 (0) 1270 580008 Fax: +44 (0) 1270 252310

Attn: Laura

Bill To: LABTECH HONG KONG LIMITED

No. 6 Anqing Street, Area B, Airport Industrial Zone, Shunyi District, Beijing 101300 P. R. China

Ship To: LABTECH HONG KONG LIMITED

No. 6 Anqing Street, Area B, Airport Industrial Zone, Shunyi District, Beijing 101300 P. R. China

Ship via: Airfreight

Quotation NO.18778

Payment: Net 30

Shipping Mark: 25HNZY051
CHANGSHA, CHINA

Shipment Expected: ASAP

Item	PN	Description	Qty.	U/M	Amount(GBP)
1	20802/4	Integra2 Integrated EA-IRMS, for isotopic and elemental analysis of 13C and 15N, or 34S. Includes standard source, high sensitivity source is available as an optional upgrade. Includes 70201-2 - analyser spares kit - filament, gaskets, pump oil, source tool kit, snoop and ferrules Includes 70501-3 - sample preparation kit - micro spatula, mirror, forceps, sample tray and syringe Includes 70401-1 - Integra2 start up consumables kit Windows based PC	1	Ea	65,000.00
2	70601-3	Sulphur option for Integra Without Dangerous goods	1	Ea	2,822.00
3	70801-4	Oxygen Upgrade Kit for Integra2- Low Temp	1	Ea	880.00
4	CON001	Customer spec pc Intel Core i5, Windows 10 ,32GB RAM, 512GB SSD 23" Monitor	1	Ea	500.00
5	WARR2	2 year warranty	1	Ea	0.00
					Total Value: GBP 69,202.00

1. Please advise us by fax or email your proposed shipment date before dispatch.

Note: Please confirm transportation before shipment.

For and on behalf of LABTECH HONG KONG LIMITED

Elena Fan

Authorized Signature

3 - MACHINE SPECIFICATION

Client:	Central South Forestry University C/O: LabTech, China		
Quotation Number:	18754	Date:	07/03/2025
Job Number:	IN1362	Date:	16/06/2025

Ref No.	Product Description	Required
HS2022 (continuous flow)		
20402/1	HS2022 Ref Gas + P Mag	
20402/2	HS2022 Ref Gas + E-Mag	
20402/3	HS2022 Ref Gas + Hd Collector + E-Mag	
20402/4	HS2022 Ref Gas + Hd Collector+E-Mag+Trap Loop	
20501	Hd Collector Assy	***
50101/2	Trapping Loop Assy	
70201-1	70201-1 Analyser spares kit – Included	
EA Units		
10202/11	iso-EArth	
10202/13	iso-EArth+	
10202/14	iso-Earth HT Furnace Upgrade	
10202/15	Iso-Earth (Without Autosampler)	
11802/1	Earth+PC	
11802/2	EArth+Clarity+PC	
70501-1	Consumables for C and N for Combustion Only	
70501-11	Simultaneous CNS Consumables Kit (approx. 10000 analysis)	
70501-3	Sample Prep Kit (Inc.)	
70501-4	GSL Combustion and Gas C & N Start Up Kit	
70501-5	Simultaneous NCS Consumables Kit	
70501-6	C&N Combustion Only Consumables Kit (approx. 5000 samples)	
70501-8	Simultaneous CNS Consumables Kit (approx. 5000 samples)	
70601-1	Sulphur Start Up Kit (Inc. with Sercon EA Systems)	
70601-2	Sulphur Consumables Kit for 1000 Samples	
70601-4	Sulphur Consumables Kit (approx. 5000 Samples)	
70801-2	Pyrolysis Consumables Kit-Low temp for Solids	
71701-1	Sub Microgram Nitrogen Analysis Start Up Kit	
SI2040	PN150 Pneumatic Autosampler with one no.50 pos, carousel	
SI2345	ME Micro Balance	
SI2737/A	Sercon 50 Position Autosampler Assy	
SI2789	Extra 50 Position Carousel	
SI2744	Automatic Liquid Sampler 15 Position	
SI2743	Liquid Sampler Mounting Plate	
Headspace Equilibration		
10202/10	XYZ Gas Sampling Control Box	
10302/1	Aqua System	
11103	µCarbs for CF	
11103-1	µCarbs for DI	
11601	200 Place Dual Heated Block 240 Volt	
11801	Sercon XYZ Gas Autosampler	
71107-1	CF oxygen-18 equilibration kit	
71107-2	DI Oxygen-18 Equilibration kit	
71201-1	CF deuterium equilibration kit	
71201-2	DI Deuterium Equilibration kit	
71801-1	WES-DI start up kit (Inc. with WES-DI)	
GC Units		
10402/2	Orchid module 1500C (Will run Hd 70901-1 C Spares Kit Inc)	
10901	Nitrogen Option (71301-1 N Spares kit Included)	
HP7890/1	Agilent 7890B GC System	
HP7890/2	Agilent 7890B GC with Auto Injector	
HP7890/3	Agilent 7890B GC + Auto Injector + 150 Vial Sample Tray	
SI2561	Nexis GC 2010 Plus	

Ref No.	Product Description	Required
---------	---------------------	----------

S027_v11

Page 1 of 6

23/10/2017

*** When ordered alone IN no. / Blue book NOT required

3 - MACHINE SPECIFICATION

Client:	Central South Forestry University C/O: LabTech, China		
Quotation Number:	18754	Date:	07/03/2025
Job Number:	IN1362	Date:	16/06/2025

GC Units *continued*

SI2564	Nexis GC-2030	
70901-3	GC Pyrolysis H and O Start Up Kit	
70901-4	GC Pyrolysis H Spares Kit	
70901-5	GC Combustion C Spares Kit (approx. 5000 samples)	
70901-6	GC Pyrolysis h Start Up Kit (approx. 5000 samples)	
71301-1	GC N Spares kit	

ABCA2 Units

10502/1	ABCA2 Sercon XYZ Autosampler	
10602	ABCA2 Gas Box	
20602	Complete Sercon ABCA2 system	
20603	ABCA IR Analyser	
70101-1	ABCA2 spares kit (<i>Inc. With ABCA2</i>)	
70101-2	ABCA2 Spares Kit for AB Analitica Systems	
70101-3	ABCA2 Autosampler Upgrade Kit	
SI1767	CVL505 – Barcode Reader	
SC9973	13C Urea breath test kit for Helicobacter pylori infection	
SC9974	Helicobacter Test kit for children age 3-11 (45mg)	

Integra2 Units

20802/4	Integra2	
70201-2	Integra2 analyser spares kit - Included	
70401-1	Integra-2 consumables (<i>Inc. With Integra2</i>)	

Trace Gas Systems

10802/4	CryoPrep	
10802/5	Automated LN refiller	
10802/6	CryoFlex (TOC)	
10802/7	CryoGas	
10802/8	CryoFlex (LA)	
SI2533	Silicon Carbide Vertical Furnace	
SI2731	213nm Laser Ablation System	
SI2733	isoScell 100 sampler chamber	
71001-1	Cryoprep Spares Kit	
71001-5	CryoGas Consumables Kit	

TOC Units

14001	Thermalox Environmental TOC-IRMS	
14002	Thermalox Environmental TOC-Standalone	
14003	Thermalox Environmental TOC-Standalone No Sampler	
SI2772	CO2 Detector, High sensitivity Option	
SI2773	NOx detector Option	
SI2774	Vacuum pump Option	
SI2775	Stirred Vial rack Option	
72101-1	TOC Start Up Kit	
72101-2	TOC Chem Pack	
72101-3	TOC & TN Instrument with heavy usage - one year	
72101-4	TOC only - economy holding	

GEO 20-22 (dual inlet)

10702	GEO 20-22 Water Equilibration System (D & 18O)	
21802/1	GEO 20-22	
21802/2	GEO 20-22 + E-Mag	
21802/3	GEO 20-22 + E-Mag +Hd	
21802/4	GEO 20-22 Inlet	

3 - MACHINE SPECIFICATION

Client:	Central South Forestry University C/O: LabTech, China		
Quotation Number:	18754	Date:	07/03/2025
Job Number:	IN1362	Date:	16/06/2025

70301-1	GEO start up kit - Included	
70301-3E	Comprehensive GEO 20-22 spares kit for 5 years operation	
HT-EA		
11501/4	HT-EA (With Autosampler)	
11501/5	HT-EA With Liquid Option	
11501/6	HT-EA Stand Alone	
11501/7	HT-EA (Without Autosampler)	
SI2622	CTC PAL3-LSI-L-S	
SI2624	CTC PAL3-Kit-Standalone	
70801-1	Pyrolysis Consumables Kit-High Temp for Solids	
71501-1	HT-EA Pyrolysis Start Up kit O and H	
71501-2	HT-EA Pyrolysis Kit for Water injections	
71501-3	HT-EA Sensitivity Start Up Kit for O & H in small samples	
71501-4	HT-EA Sensitivity Consumables Kit for O & H in small samples	
71501-5	HT-EA Sensitivity Kit for O&H in Small Samples (5000 samples)	
71501-7	HT-EA Kit for Room Temp Equilibration	
71501-8	HT-EA Pyrolysis Kit for Water injections (approx. 1 year kit)	

Ref No.	Product Description	Required
Upgrades		
12101	Evacuation System ***	
12102	2 Ref Gas+ Dilutor Box ***	
20101	Standard Source Assy ***	
20102	HS Source Assy ***	
20301	Standard Collector Assy ***	
30202/1	System Control Upgrade (Pmag) installation not included	Edwards Or Pfeiffer
30202/2	System Controller Upgrade (Emag) installation not included	Edwards Or Pfeiffer
10202/7	Preparation PCB Upgrade	
30202/3	Hijack Box ***	
30302/4	nEXT85 Edwards Pumping System Upgrade Single - installation not included	
30302/5	nEXT85 Edwards Pumping System Upgrade Dual - installation not included	
30302/6	nEXT85 Edwards Pumping System Upgrade Triple- installation not included	
SI2542	ABCA2 Upgrade Kit	
CALLISTOW	Callisto software upgrade ***	
Specials		
11702	LC-I interface	
12202	Leco Interface Unit	
SI2320	Agilent 1260 Series Refractive Index Detection HPLC System	
71401-1	ABCA-G Start Up/Spares Kit	
71401-2	ABCA-G Equilibration Spares Kit - Room temp O only	
71401-3	ABCA-G Equilibration Spares Kit - Room Temp O & H	
71601-1	LC Interface Spares Kit	
ADDITIONS		
PC	Windows 7 based PC System with Callisto Software	
Printer	Printer	
INST1	Installation and start up training at customer's site by Sercon Engineer	
INST2	Installation and start up training at customer's site by Agents Engineer	
WARR	1 year warranty	
WARR2	2 year warranty	
TRAIN	Training in UK for one person for 3 days	
SCONTRACT	Maintenance Service Contract – 1 year	
SPECIFICATION REQUIREMENTS		
Voltage	Volts 50/60Hzs	

3 - MACHINE SPECIFICATION

Client:	Central South Forestry University C/O: LabTech, China		
Quotation Number:	18754	Date:	07/03/2025
Job Number:	IN1362	Date:	16/06/2025

Plug Type	SI2510	3 Metre USA/Japan plug to IEC C 13	
	SI2511	2.5 Metre Australian plug to IEC C 13	
	SI2512	3 Metre UK 13 amp Fused plug to IEC C 13	
	SI2523	3metre Schuko (Euro) Plug to IEC C13 Socket	
	SI2514	2.5 metre Chinese mains to IEC C13	
Pallet	Stamped for export heat treated for 30 min @ 56 °C		
	Complies with all Export Requirements		

4 - LOG SHEET

Client:			
Quotation Number:		Date:	
Job Number:		Date:	

Date	Action Taken	By Whom
12/3/25	UNIT P.A.T TESTED.	OLI
17/3/25	UNIT PUMPED DOWN - WITHOUT 6 THOU PASSUP 5 E-8. BUT THEN CONTROLLER WENT OFF	
18/3/25	UNIT PUMPED DOWN WITH 6 THOU PASSUP 5 E-8. BUT THEN CONTROLLER WENT OFF.	
26.3.25	System turned on and pump out down. All crimps set and pressure tested Reaction tubes and GC left heated overnight	AD/TI
27.3.25	Isolation valve crimped. Ran Samples for tuning. left a precision of 400/100 overnight	AD/TI
28.3.25	Precision passed, ran linearity, failed. retuned and ran a new precision. CO ₂ shifted back, baking GC @ 140°C	AD
31.3.25	Blanks high, leak at water trap, swapped red caps for new supply. Pressure test passed Running a linearity (prec ran as usual) and passed	AD
1.4.25	linearity passed	
2.4.25	Ran Cystine against itself and against GMA, Very close to correct delta's. running again with half the amount of GMA. Blank high, but leak check passes, checking tubes	AD

Date	Action Taken	By Whom
------	--------------	---------



sercon
innovators in isotopes

Sercon Ltd, Europa House, Coppicemere Drive,
Crewe, Cheshire CW1 6GZ, UK
Tel: +44 (0)1270 580008
Company Registration No. 4290072 Vat No. 771 1100 74

Test acceptance - 6

6 - TEST / ACCEPTANCE SPECIFICATION

Customer:			
Quotation Number:		Date:	
Job Number:		Date:	

HS20-22			
Reference Gas Precision (n=5) x5 Ref gas pulses of 5 seconds to be used.			
¹³ C	(CO ₂)	0.06	
¹⁵ N		0.06	
¹⁸ O	(CO ₂)	0.06	
³⁴ S	(SO ₂)	0.1	
² H		0.4	

Reference Gas Linearity ref gas pulses ranging from 1-second trough to 10 seconds			
¹³ C	≤0.02‰/nA	STD Deviation <0.30	
¹⁵ N	≤0.02‰/nA	STD Deviation <0.30	

IsoEarth & PyroEarth					
Combustion (n=5) (For NCS use 0.3mg of Cystine) (Pyrolysis use silver caps) 0.4mg of Cystine for ¹⁸ O and 4 mg of Cystine for H ₂					
¹³ C	100 µg C	0.1		10 µg C	0.8
¹⁵ N	100 µg N	0.15		10 µg N	1.5
³⁴ S	25 µg S	0.3		10 µg S	1.0
Dual	400 µg C	0.1		100 µg N	0.15
¹⁸ O	100 µg or 1 µl water	0.3			
² H	200 µg H or 1 µl water	1.0			

6 - TEST / ACCEPTANCE SPECIFICATION

Customer:			
Quotation Number:		Date:	
Job Number:		Date:	

EArth			
Solid specs to be completed using 2mg methionine calibration against 1,2,3mg (n=3) cystine. Combustion (n=5)			
		Accuracy (avg)	Precision (std dev)
¹³ C		+/- 0.4%	0.2
¹⁵ N		+/- 0.4%	0.2
² H		+/- 0.4%	0.2
³⁴ S		+/- 0.4%	0.2
¹⁸ O	Pyrolysis	+/- 0.4%	0.2
Liquid injections n=5 using 5µl methanol calibrated against 2,4,6µl injections of ethanol			
		Accuracy (avg)	Precision (std dev)
¹⁸ O	Pyrolysis	+/- 0.4%	0.2

Integra2			
Specs to be completed using AS beet mix 400C/100N for precision & 250C/100N for linearity			
¹³ C	100 µg C	0.2	
¹⁵ N	100 µg N	0.3	
Dual	100 µg N & 400 µg C	0.4 & 0.15	CO ₂ = 0.076 N ₂ = 0.067
³⁴ S	30 µg S	0.5	
¹⁸ O	100 µg O	0.5	
Dual	100 µg N & 250 µg C	0.5	CO ₂ = 0.43 N ₂ = 0.09
HT-EA (Pyrolysis use silver caps) 0.4mg of Cystine for ¹⁸ O and 4 mg of Cystine for H ₂			
Pyrolysis (n=5)			
¹⁸ O	100 µg O or 1 µl water	0.3	
² H	200 µg H or 1 µl water	1.0	

6 - TEST / ACCEPTANCE SPECIFICATION

Customer:			
Quotation Number:		Date:	
Job Number:		Date:	

Aqua, μ Carbs & GasBox					
Equilibrations (n=5)					
^{18}O	0.2 ml water	0.1			
^2H	0.2 ml water	2.0			
Carbonates (n=5)					
^{13}C	100 μg CaCO_3	0.08			
^{18}O	100 μg CaCO_3	0.1			
Gas Analysis (n=5)					
^{13}C	5% CO_2	0.1		Air CO_2	0.1
^{18}O	5% CO_2	0.1		Air CO_2	0.2
^2H	0.5 ml H_2	1.5			
^{15}N	0.5 ml N_2	0.1		100 ppm N_2O	0.5
^{13}C	Breath (12 ml)	0.1			
DIC (n=5)					
^{13}C	15ppm (0.5ml)	0.15			

LC-I					
(n=5)					
^{13}C		0.2			

GC Orchid					
Combustion (n=3 for each analyte)					
^{13}C	3 Alkanes	100 ng C	0.2	10 ng C	0.5
^{15}N	2 Amines	100 ng N	0.5	14 ng N	1.5

6 - TEST / ACCEPTANCE SPECIFICATION

Customer:			
Quotation Number:		Date:	
Job Number:		Date:	

Pyrolysis (n=3 for each analyte)						
² H	3 Alkanes	250 ng H	3.0		25 ng H	10.0
¹⁸ O	3 FAMES	100 ng O	0.5		20 ng O	2.0
Water pyrolysis (n=5)						
² H	0.1 µl water	3.0				
¹⁸ O	0.1 µl water	0.4				

CryoRange						
(n=5)						
¹³ C	Air CO ₂ (12ml)	0.1		Air CH ₄ (125ml)	0.2	
¹⁵ N	Air N ₂ (1ml)	0.1		Air N ₂ O (125ml)	0.4	
¹⁸ O	Air CO ₂ (12ml)	0.15		Air N ₂ O (125ml)	0.7	
¹⁸ O	Air O ₂ (1ml)	0.1				
² H	Air CH ₄	3.0				
¹³ C	TOC (5ppm)	0.4				
¹³ C	Laser					

Geo2			
Internal precision (no change overs)			
² H	200 Bar µl	0.15	
¹³ C	100 Bar µl	0.01	
¹⁵ N	100 Bar µl	0.01	
¹⁸ O	100 Bar µl	0.015	
³⁴ S	100 Bar µl	0.015	
Single cold finger volumes (10 change overs)			

6 - TEST / ACCEPTANCE SPECIFICATION

Customer:			
Quotation Number:		Date:	
Job Number:		Date:	

¹³ C	20 Bar µl	0.01	
¹⁸ O	20 Bar µl	0.015	
¹³ C	3 Bar µl	0.02	
¹⁸ O	3 Bar µl	0.03	
CAPS (n=5)			
¹³ C	50 µg CaCO ₃	0.03	
¹⁸ O	50 µg CaCO ₃	0.06	
WES (n=5)			
² H		1.0	
¹⁸ O		0.05	
Manifold (10 tubes)			
² H		0.5	
¹³ C		0.02	
¹⁵ N		0.02	
¹⁸ O		0.03	
³⁴ S		0.03	

Special tests



Sercon Ltd, Europa House, Coppicemere Drive,
Crewe, Cheshire CW1 6GZ, UK
Tel: +44 (0)1270 580008
Company Registration No. 4290072 Vat No. 771 1100 74

Test acceptance - 6

6 - TEST / ACCEPTANCE SPECIFICATION

Customer:			
Quotation Number:		Date:	
Job Number:		Date:	

Sales Signature: 	Date: 3/07/25
Test Signature: 	Date: 3/07/25

7: ASSEMBLY AND TEST RECORD

Client:			
Quotation Number:		Date:	
Job Number:		Date:	

PRODUCT:

INTEGRA 2

SERIAL NUMBER:

20802/4

 ASSEMBLED:
 (Name)

Oli

[Signature]

Date:

12/3/25

 PRE-TEST
 INSPECTION:
 (Production name)

P. Mayo

Date:

18/3/25

 PRE-TEST
 INSPECTION:
 (Test name)

T. Iordanidis

Date:

18/3/25

 SYSTEM TEST:
 (Name)

[Signature]

Date:

8/5/25

 PRE-DESPATCH
 INSPECTION:
 (Name)

[Signature]

Date:

2/6/25

P.A.T. Test results

 In compliance with Low Voltage directive 73/23/EEC
 To Low Voltage standard EN61010 – 1:1993

Earth Continuity

0.072 Ω

Insulation

7300 $M\Omega$

Earth Leakage Current

1.90mA mA. @ 0.12 kVA.

Tested by:

(Oli)

Date:

12/3/25

8 - Serial Numbers Record Document

Client:	Central South Forestry University C/O: LabTech, China		
Quotation Number:	18754	Date:	07/03/2025
Job Number:	IN1362	Date:	16/06/2025

Units	S/N
MS/ABCA	
Integra 2	20802/4/2025/1362
Prep Unit 1	
Prep Unit 2/Misc	
Electronics Upgrade	
CETAC	
Gas Box	
Autosampler	42694
BCR	
Dual Heater Block	

Pump Kit and Optics	S/N
Turbo pump S	250071198
Turbo pump A	
Rotary Pump	250065857
Voltage	110 / 240
Gauge S	250009658
Gauge A	
Source	20101-251
Collector 1	20301-167
Collector 2	
Magnet	18/7/23

Parts	S/N
Fan PSU	
24V PSU	

PCBs	S/N
Electronics Box	20190006
Prep PCB	
VFC PCB	
Watchdog PCB	
Head Amp PCB 1	
Head Amp PCB 2	24025
Electromagnet RP	

Furnaces	S/N
Combustion	
Reduction	
HT-EA	
GC Oven	
Furnace 1	
Furnace 2	

PC	S/N
PC	B2W7T784
Monitor	H924Q44
Printer	

Key

Pumps:	E: Edwards
	P: Pfeiffer
Furnaces:	V/S: Vecstar
	E: Elite
Magnet:	E: Electromagnet
	S: Standard
Turbos & Gauges:	S: Source
	A: Analyser
Collectors:	Std: Standard
	S: Sulphur
	HD: Hydrogen Deuterium

Additional System Info:



Packing List

Customer Reference: 20250500013

Shipper / Exporter:

Sercon Limited
Europa House
Coppicemere Drive
Crewe
CW1 6GZ
UK

Tel: +44 (0)1270 580008
Reference: IN1362

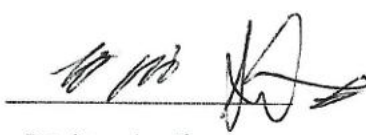
Consignee:

LabTech Corporation Inc
No. 6 Anqing Street, Area B
Airport Industrial Park
Shunyi District
Beijing, 101300
P.R. CHINA

Contact Name: Elena Fan
Contact Tel Number: 0086-10-8048 6450

BOX: 1

Quantity	Product Code	Product Description	Serial Number: if available
☒ 1	20802/4	Integra	20802/4/2025/1362

Packed by: 

Date: 03/07/25

Claims for damage or shortages must be made within 3 days



S035b_v6

Sercon Limited
Europa House, Coppicemere Drive
Crewe, Cheshire
CW1 6GZ, U.K.
Company Registration No.: 4290072

Phone: +44 (0) 1270 580008

sales@sercongroup.com
www.sercongroup.com
VAT No.: 771 1100 74
19/01/2017

Packing List

Customer Reference: 20250500013

Shipper / Exporter:

Sercon Limited
Europa House
Coppicemere Drive
Crewe
CW1 6GZ
UK

Tel: +44 (0)1270 580008
Reference: IN1362

Consignee:

LabTech Corporation Inc
No. 6 Anqing Street, Area B
Airport Industrial Park
Shunyi District
Beijing, 101300
P.R. CHINA

Contact Name: Elena Fan
Contact Tel Number: 0086-10-8048 6450

BOX: 2

Quantity	Product Code	Product Description	Serial Number: if available
☒ 1	70501-3	Sample prep kit	
☒ 1	70401-1	Integra 2 startup kit	
☒ 1	70201-2	Integra 2 analyser spares kit	
☒ 1	SI0101	Magnet	18/7/23
☒ 1	SI0103	Autosampler	42694
☒ 1	SI0103	Rotary pump	2500658557
☒ 1	SI2391	Monitor	H924G44
☒ 1	CON 002	Dell 7020 small form factor	B2W7T74
☒ 1	70601-3	Sulphur option for integra	
☒ 1	70801-4	Oxygen upgrade	
☒ 1		Software backup USB (inside PC box)	
☒ 1		Windows recovery USB (inside PC box)	

Packed by: 

Date: 03/07/25

Claims for damage or shortages must be made within 3 days