

CONTENTS PAGE

Client:	Weifang Nongfeng Testing Co., Ltd C/O: LabTech, China		
Quotation Number:	16668	Date:	17/10/2022
Job Number:	IN1294	Date:	25/10/2022

- ☒ 1 QUOTATION
- ☒ 2 CUSTOMER ORDER
- ☒ 3 MACHINE SPECIFICATION DOCUMENT_S027_v11
- ☒ 4 LOG SHEET_S028_v3
- ☒ 5 FAULT REPORTING DOCUMENT_S029_v3
- ☒ 6 TEST ACCEPTANCE SPECIFICATION_S030_v9
- ☒ 7 ASSEMBLY AND TEST RECORD_S023_v5
- ☒ 8 SERIAL NUMBER FORM_S053_v4
- ☒ 9 TEST DATA
- ☒ 10 PACKING LIST_S035b_v6
- ☐ 11 INSTALLATION REPORT_S057_v6
- ☐ 12 INSTALLATION ACCEPTANCE AND WARRANTY CERTIFICATE_S071_v3

ADDITIONAL:

Payment Terms:
 50% / 40% / 10%
Other: please specify _____

Penalty Clause Yes / No

Please specify: _____

Blue Book Completion:

Signed	<i>D. O'Connor</i>	Date	28/2/23
Verified	<i>S. Mahon</i>	Date	28/2/23

3 - MACHINE SPECIFICATION

Client:	Weifang Nongfeng Testing Co., Ltd C/O: LabTech, China		
Quotation Number:	16668	Date:	17/10/2022
Job Number:	IN1294	Date:	25/10/2022

Ref No.	Product Description	Required
HS2022 (continuous flow)		
20402/1	HS2022 Ref Gas + P Mag	
20402/2	HS2022 Ref Gas + E-Mag	
20402/3	HS2022 Ref Gas + Hd Collector + E-Mag	
20402/4	HS2022 Ref Gas + Hd Collector+E-Mag+Trap Loop	
20501	Hd Collector Assy	***
50101/2	Trapping Loop Assy	
70201-1	70201-1 Analyser spares kit - Included	
EA Units		
10202/11	iso-EArth	
10202/13	iso-EArth+	
10202/14	iso-Earth HT Furnace Upgrade	
10202/15	Iso-Earth (Without Autosampler)	
11802/1	Earth+PC	
11802/2	EArth+Clarity+PC	
70501-1	Consumables for C and N for Combustion Only	
70501-11	Simultaneous CNS Consumables Kit (approx. 10000 analysis)	
70501-3	Sample Prep Kit (Inc.)	
70501-4	GSL Combustion and Gas C & N Start Up Kit	
70501-5	Simultaneous NCS Consumables Kit	
70501-6	C&N Combustion Only Consumables Kit (approx. 5000 samples)	
70501-8	Simultaneous CNS Consumables Kit (approx. 5000 samples)	
70601-1	Sulphur Start Up Kit (Inc. with Sercon EA Systems)	
70601-2	Sulphur Consumables Kit for 1000 Samples	
70601-4	Sulphur Consumables Kit (approx. 5000 Samples)	
70801-2	Pyrolysis Consumables Kit-Low temp for Solids	
71701-1	Sub Microgram Nitrogen Analysis Start Up Kit	
SI2040	PN150 Pneumatic Autosampler with one no.50 pos, carousel	
SI2345	ME Micro Balance	
SI2737/A	Sercon 50 Position Autosampler Assy	
SI2789	Extra 50 Position Carousel	
SI2744	Automatic Liquid Sampler 15 Position	
SI2743	Liquid Sampler Mounting Plate	
Headspace Equilibration		
10202/10	XYZ Gas Sampling Control Box	
10302/1	Aqua System	
11103	µCarbs for CF	
11103-1	µCarbs for DI	
11601	200 Place Dual Heated Block 240 Volt	
11801	Sercon XYZ Gas Autosampler	
71107-1	CF oxygen-18 equilibration kit	
71107-2	DI Oxygen-18 Equilibration kit	
71201-1	CF deuterium equilibration kit	
71201-2	DI Deuterium Equilibration kit	
71801-1	WES-DI start up kit (Inc. with WES-DI)	
GC Units		
10402/2	Orchid module 1500C (Will run Hd 70901-1 C Spares Kit Inc)	
10901	Nitrogen Option (71301-1 N Spares kit Included)	
HP7890/1	Agilent 7890B GC System	
HP7890/2	Agilent 7890B GC with Auto Injector	
HP7890/3	Agilent 7890B GC + Auto Injector + 150 Vial Sample Tray	
SI2561	Nexis GC 2010 Plus	

3 - MACHINE SPECIFICATION

Client:	Weifang Nongfeng Testing Co., Ltd C/O: LabTech, China		
Quotation Number:	16668	Date:	17/10/2022
Job Number:	IN1294	Date:	25/10/2022

Ref No.	Product Description	Required
GC Units <i>continued</i>		
SI2564	Nexis GC-2030	
70901-3	GC Pyrolysis H and O Start Up Kit	
70901-4	GC Pyrolysis H Spares Kit	
70901-5	GC Combustion C Spares Kit (approx. 5000 samples)	
70901-6	GC Pyrolysis h Start Up Kit (approx. 5000 samples)	
71301-1	GC N Spares kit	
ABCA2 Units		
10502/1	ABCA2 Sercon XYZ Autosampler	
10602	ABCA2 Gas Box	
20602	Complete Sercon ABCA2 system	
20603	ABCA IR Analyser	
70101-1	ABCA2 spares kit (Inc. With ABCA2)	
70101-2	ABCA2 Spares Kit for AB Analitica Systems	
70101-3	ABCA2 Autosampler Upgrade Kit	
SI1767	CVL505 - Barcode Reader	
SC9973	13C Urea breath test kit for Helicobacter pylori infection	
SC9974	Helicobacter Test kit for children age 3-11 (45mg)	
Integra2 Units		
20802/4	Integra2	
70201-2	Integra2 analyser spares kit - Included	
70401-1	Integra-2 consumables (Inc. With Integra2)	
Trace Gas Systems		
10802/4	CryoPrep	
10802/5	Automated LN refiller	
10802/6	CryoFlex (TOC)	
10802/7	CryoGas	
10802/8	CryoFlex (LA)	
SI2533	Silicon Carbide Vertical Furnace	
SI2731	213nm Laser Ablation System	
SI2733	isoScell 100 sampler chamber	
71001-1	Cryoprep Spares Kit	
71001-5	CryoGas Consumables Kit	
TOC Units		
14001	Thermalox Environmental TOC-IRMS	
14002	Thermalox Environmental TOC-Standalone	
14003	Thermalox Environmental TOC-Standalone No Sampler	
SI2772	CO2 Detector, High sensitivity Option	
SI2773	NOx detector Option	
SI2774	Vacuum pump Option	
SI2775	Stirred Vial rack Option	
72101-1	TOC Start Up Kit	
72101-2	TOC Chem Pack	
72101-3	TOC & TN Instrument with heavy usage - one year	
72101-4	TOC only - economy holding	
GEO 20-22 (dual inlet)		
10702	GEO 20-22 Water Equilibration System (D & 18O)	
21802/1	GEO 20-22	
21802/2	GEO 20-22 + E-Mag	

3 - MACHINE SPECIFICATION

Client:	Weifang Nongfeng Testing Co., Ltd C/O: LabTech, China		
Quotation Number:	16668	Date:	17/10/2022
Job Number:	IN1294	Date:	25/10/2022

21802/3	GEO 20-22 + E-Mag +Hd	
21802/4	GEO 20-22 Inlet	
70301-1	GEO start up kit - Included	
70301-3E	Comprehensive GEO 20-22 spares kit for 5 years operation	
HT-EA		
11501/4	HT-EA (With Autosampler)	
11501/5	HT-EA With Liquid Option	
11501/6	HT-EA Stand Alone	
11501/7	HT-EA (Without Autosampler)	
SI2622	CTC PAL3-LSI-L-S	
SI2624	CTC PAL3-Kit-Standalone	
70801-1	Pyrolysis Consumables Kit-High Temp for Solids	
71501-1	HT-EA Pyrolysis Start Up kit O and H	
71501-2	HT-EA Pyrolysis Kit for Water injections	
71501-3	HT-EA Sensitivity Start Up Kit for O & H in small samples	
71501-4	HT-EA Sensitivity Consumables Kit for O & H in small samples	
71501-5	HT-EA Sensitivity Kit for O&H in Small Samples (5000 samples)	
71501-7	HT-EA Kit for Room Temp Equilibration	
71501-8	HT-EA Pyrolysis Kit for Water injections (approx. 1 year kit)	

Ref No.	Product Description	Required
Upgrades		
12101	Evacuation System ***	
12102	2 Ref Gas+ Dilutor Box ***	
20101	Standard Source Assy ***	
20102	HS Source Assy ***	
20301	Standard Collector Assy ***	
30202/1	System Control Upgrade (Pmag) installation not included	Edwards Or Pfeiffer
30202/2	System Controller Upgrade (Emag) installation not included	Edwards Or Pfeiffer
10202/7	Preparation PCB Upgrade	
30202/3	Hijack Box ***	
30302/4	nEXT85 Edwards Pumping System Upgrade Single - installation not included	
30302/5	nEXT85 Edwards Pumping System Upgrade Dual - installation not included	
30302/6	nEXT85 Edwards Pumping System Upgrade Triple- installation not included	
SI2542	ABCA2 Upgrade Kit	
CALLISTOW	Callisto software upgrade ***	
Specials		
11702	LC-I interface	
12202	Leco Interface Unit	
SI2320	Agilent 1260 Series Refractive Index Detection HPLC System	
71401-1	ABCA-G Start Up/Spares Kit	
71401-2	ABCA-G Equilibration Spares Kit - Room temp O only	
71401-3	ABCA-G Equilibration Spares Kit - Room Temp O & H	
71601-1	LC Interface Spares Kit	
ADDITIONS		
PC	Windows 7 based PC System with Callisto Software	✓
Printer	Printer	
INST1	Installation and start up training at customer's site by Sercon Engineer	
INST2	Installation and start up training at customer's site by Agents Engineer	✓
WARR	1 year warranty	
WARR2	2 year warranty	✓
TRAIN	Training in UK for one person for 3 days	
SCONTRACT	Maintenance Service Contract – 1 year	

3 - MACHINE SPECIFICATION

Client:	Weifang Nongfeng Testing Co., Ltd C/O: LabTech, China		
Quotation Number:	16668	Date:	17/10/2022
Job Number:	IN1294	Date:	25/10/2022

SPECIFICATION REQUIREMENTS			
Voltage	Volts 50/60Hzs	220V	
Plug Type	SI2510	3 Metre USA/Japan plug to IEC C 13	✓
	SI2511	2.5 Metre Australian plug to IEC C 13	
	SI2512	3 Metre UK 13 amp Fused plug to IEC C 13	
	SI2523	3metre Schuko (Euro) Plug to IEC C13 Socket	
	SI2514	2.5 metre Chinese mains to IEC C13	
Pallet	Stamped for export heat treated for 30 min @ 56 °C		✓
	Complies with all Export Requirements		✓

4 - LOG SHEET

Client:	Weifang Nongfeng Testing Co., Ltd C/O: LabTech, China		
Quotation Number:	16668	Date:	17/10/2022
Job Number:	IN1294	Date:	25/10/2022

Date	Action Taken	By Whom
3/1/22	Instrument put into pump down cycle Combustion & reduction tubes installed CO₂ & Emg. satb trap installed System hooked up to main H₂O₂ line And system pressure tested. Large leak on auto-sampler (S) Greik screw was not sealed and bolts not torqued down. This was rectified until the A.S. became leak tight. Struggled to get sampler holes to line up? Resumed pressure testing. Which on the bottom both furnaces were too short and misaligned with carrier plate. Broke the reduction tube trying to fit the tubes in situ. Pressure tested O.K. Cimps installed and flows set. Flow meter calibrated G.C baked out 140°C overnight G.C reduced to 65°C + furnaces set to working temps. Filament conditioned	SB.
4/1/23	G.C baked out 140°C overnight G.C reduced to 65°C + furnaces set to working temps. Filament conditioned	SB.
N.B	Flight tube had to be baked out due to contamination. Vacuum degraded constantly after filament conditioning. Isolation valve stopped closing. Valve removed and spring stretched & cap hole enlarged. Source flashing over @ 3Kv. Trying to condition source to prevent	
5/1/23	Vacuum degrades with increases in H.T. trying to condition the source. Source removed and a piece of wire fell from the source. The rest of the source was in good shape	SB

4 - LOG SHEET


Client:	Weifang Nongfeng Testing Co., Ltd C/O: LabTech, China		
Quotation Number:	16668	Date:	17/10/2022
Job Number:	IN1294	Date:	25/10/2022

Date	Action Taken	By Whom
5/1/23	and the filament was in good condition Source reinstalled and system put into pump down cycle. Tuned up source to Garching parameters just a small fluctuation in vacuum level noted. Increased H.T to 3.5KV Vacuum level remained stable. Tuned up for CO ₂	SB
	Software when first opened has Timeout message near beams. It disappears when a value is operated. Wayne to investigate Monday	SM
	New ion source fitted. Instrument pumping down as of 11am.	
	Filament conditioned fine	SM.
9-1-23	Excessive emission on Standby setting. Source removed & checked filament slightly out of line in one magnet too far back in the holder. Source reinstalled and system pumped down.	SB
16-1-23	Source rebuilt and reinstalled. Pumped down and conditioned. Source still arcing over on N ₂ settings Changed source cable and arcing stopped. Source now only arcs when isolation valve open? System leak checked O.K. Isolation valve changed due to poor performance. 4.5e11. Bg signal is closed. 4.8e11 " " " open.	SB




Client:	Weifang Nongfeng Testing Co., Ltd C/O: LabTech, China		
Quotation Number:	16668	Date:	17/10/2022
Job Number:	IN1294	Date:	25/10/2022

S028_v3

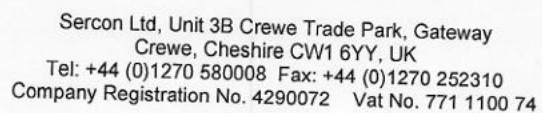

Client:	Weifang Nongfeng Testing Co., Ltd C/O: LabTech, China		
Quotation Number:	16668	Date:	17/10/2022
Job Number:	IN1294	Date:	25/10/2022

S028_v3

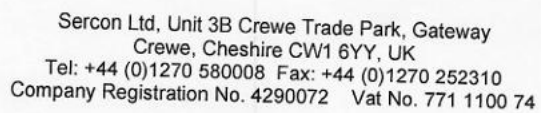



Test acceptance - 6

Sercon Ltd, Unit 3B Crewe Trade Park, Gateway
Crewe, Cheshire CW1 6JT, UK
Tel: +44 (0)1270 580008 Fax: +44 (0)1270 252310
Company Registration No. 4290072 Vat No. 771 1100 74

6 - TEST / ACCEPTANCE SPECIFICATION

Customer:	Weifang Nongfeng Testing Co., Ltd C/O: LabTech, China		
Quotation Number:	16668	Date:	17/10/2022
Job Number:	IN1294	Date:	25/10/2022

Special tests

Sales Signature: <u>Joe Connor</u>	Date: <u>1/2/23</u>
Test Signature: <u>S. Masely</u>	Date: <u>1/2/23</u>

6 - TEST / ACCEPTANCE SPECIFICATION

Customer:	Weifang Nongfeng Testing Co., Ltd C/O: LabTech, China		
Quotation Number:	16668	Date:	17/10/2022
Job Number:	IN1294	Date:	25/10/2022

HS20-22

Reference Gas Precision (n=5)

¹³ C	(CO ₂)	0.06	
¹⁵ N		0.06	
¹⁸ O	(CO ₂)	0.06	
³⁴ S	(SO ₂)	0.1	
² H		0.4	

Reference Gas Linearity

¹³ C	≤0.02‰/nA	
¹⁵ N	≤0.02‰/nA	

IsoEarth & IsoEarth+

Combustion (n=5)

¹³ C	100 µg C	0.1		10 µg C	0.8	
¹⁵ N	100 µg N	0.15		10 µg N	1.5	
³⁴ S	25 µg S	0.3		10 µg S	1.0	
Dual	400 µg C	0.1		100 µg N	0.15	

Integra2

¹³ C	100 µg C	0.2	0.03	
¹⁵ N	100 µg N	0.3	0.07 0.07	
Dual	100 µg N & 400 µg C	0.4 & 0.15	0.18 & 0.05	
³⁴ S	30 µg S	0.5		
¹⁸ O	100 µg O	0.5		

Lin³
15N

0.37
0.26

S030_v9

07/10/2019

6 - TEST / ACCEPTANCE SPECIFICATION

Customer:	Weifang Nongfeng Testing Co., Ltd C/O: LabTech, China		
Quotation Number:	16668	Date:	17/10/2022
Job Number:	IN1294	Date:	25/10/2022

HT-EA			
Pyrolysis (n=5)			
¹⁸ O	100 µg O or 1 µl water	0.3	
² H	200 µg H or 1 µl water	1.0	

Aqua, µCarbs & GasBox					
Equilibrations (n=5)					
¹⁸ O	0.2 ml water	0.1			
² H	0.2 ml water	2.0			
Carbonates (n=5)					
¹³ C	100 µg CaCO ₃	0.08			
¹⁸ O	100 µg CaCO ₃	0.1			
Gas Analysis (n=5)					
¹³ C	5% CO ₂	0.1		Air CO ₂	0.1
¹⁸ O	5% CO ₂	0.1		Air CO ₂	0.2
² H	0.5 ml H ₂	1.5			
¹⁵ N	0.5 ml N ₂	0.1		100 ppm N ₂ O	0.5
¹³ C	Breath (12 ml)	0.1			
DIC (n=5)					
¹³ C	15ppm (0.5ml)	0.15			

LC-I			
(n=5)			
¹³ C		0.2	

6 - TEST / ACCEPTANCE SPECIFICATION

Customer:	Weifang Nongfeng Testing Co., Ltd C/O: LabTech, China		
Quotation Number:	16668	Date:	17/10/2022
Job Number:	IN1294	Date:	25/10/2022

GC Orchid						
Combustion (n=3 for each analyte)						
¹³ C	3 Alkanes	100 ng C	0.2		10 ng C	0.5
¹⁵ N	2 Amines	100 ng N	0.5		14 ng N	1.5
Pyrolysis (n=3 for each analyte)						
² H	3 Alkanes	250 ng H	3.0		25 ng H	10.0
¹⁸ O	3 FAMES	100 ng O	0.5		20 ng O	2.0
Water pyrolysis (n=5)						
² H	0.1 µl water	3.0				
¹⁸ O	0.1 µl water	0.4				

CryoRange						
(n=5)						
¹³ C	Air CO ₂ (12ml)	0.1		Air CH ₄ (125ml)	0.2	
¹⁵ N	Air N ₂ (1ml)	0.1		Air N ₂ O (125ml)	0.4	
¹⁸ O	Air CO ₂ (12ml)	0.15		Air N ₂ O (125ml)	0.7	
¹⁸ O	Air O ₂ (1ml)	0.1				
² H	Air CH ₄	3.0				
¹³ C	TOC (5ppm)	0.4				
¹³ C	Laser					

6 - TEST / ACCEPTANCE SPECIFICATION

Customer:	Weifang Nongfeng Testing Co., Ltd C/O: LabTech, China		
Quotation Number:	16668	Date:	17/10/2022
Job Number:	IN1294	Date:	25/10/2022

Geo2			
Internal precision (no change overs)			
² H	200 Bar µl	0.15	
¹³ C	100 Bar µl	0.01	
¹⁵ N	100 Bar µl	0.01	
¹⁸ O	100 Bar µl	0.015	
³⁴ S	100 Bar µl	0.015	
Single cold finger volumes (10 change overs)			
¹³ C	20 Bar µl	0.01	
¹⁸ O	20 Bar µl	0.015	
¹³ C	3 Bar µl	0.02	
¹⁸ O	3 Bar µl	0.03	
CAPS (n=5)			
¹³ C	50 µg CaCO ₃	0.03	
¹⁸ O	50 µg CaCO ₃	0.06	
WES (n=5)			
² H		1.0	
¹⁸ O		0.05	
Manifold (10 tubes)			
² H		0.5	
¹³ C		0.02	
¹⁵ N		0.02	
¹⁸ O		0.03	
³⁴ S		0.03	

7: ASSEMBLY AND TEST RECORD

Client:	Weifang Nongfeng Testing Co., Ltd C/O: LabTech, China		
Quotation Number:	16668	Date:	17/10/2022
Job Number:	IN1294	Date:	25/10/2022

PRODUCT:

INTEGRA 2

SERIAL NUMBER:

2080214-2023-1294

 ASSEMBLED:
 (Name)

Jeanne Ham

Date:

20/12/2022

 PRE-TEST
 INSPECTION:
 (Production name)

Jeanne Ham

Date:

20/12/2022

 PRE-TEST
 INSPECTION:
 (Test name)

Steven Bennett

Date:

3/1/2023

 SYSTEM TEST:
 (Name)

~~Jeanne Ham~~ Steven Bennett

Date:

25/1/2023

 PRE-DESPATCH
 INSPECTION:
 (Name)

Jeanne Ham

Date:

27/1/23

P.A.T. Test results

 In compliance with Low Voltage directive 73/23/EEC
 To Low Voltage standard EN61010 – 1:1993

Earth Continuity

0.077 Ω

Insulation

7300 $M\Omega$

Earth Leakage Current

1.15 mA. @ 0.05 kVA.

Tested by:

P. Hays

Date:

19/12/22

8 - Serial Numbers Record Document

Client:	Weifang Nongfeng Testing Co., Ltd C/O: LabTech, China		
Quotation Number:	16668	Date:	17/10/2022
Job Number:	IN1294	Date:	25/10/2022

Units	S/N
MS/ABCA	
Integra 2	20802/4-2023/1294
Prep Unit 1	
Prep Unit 2/Misc	
Electronics Upgrade	
CETAC	
Gas Box	
Autosampler	42666
BCR	
Dual Heater Block	

Pump Kit and Optics	S/N
^{AGILENT} Turbo pump S	1T22436111
Turbo pump A	
^{AGILENT} Rotary Pump	JP22200036
Voltage	110 / 240
Gauge S	L12238F087
Gauge A	
Source	20101-203
Collector 1	20301-124
Collector 2	
Magnet	21/12/22

Parts	S/N
Fan PSU	
24V PSU	

PCBs	S/N
Electronics Box	20230001
Prep PCB	
VFC PCB	
Watchdog PCB	
Head Amp PCB 1	22023
Head Amp PCB 2	
Electromagnet RP	

Furnaces	S/N
Combustion	
Reduction	
HT-EA	
GC Oven	
Furnace 1	
Furnace 2	

PC	S/N
PC	8V2DQR3
Monitor	3QMP8H3
Printer	

Key

Pumps:	E: Edwards P: Pfeiffer
Furnaces:	V/S: Vecstar E: Elite
Magnet:	E: Electromagnet S: Standard
Turbos & Gauges:	S: Source A: Analyser
Collectors:	Std: Standard S: Sulphur HD: Hydrogen Deuterium

Additional System Info:



Packing List

Customer Reference: 20221000050

Shipper / Exporter:

Sercon Limited
Unit 3B Crewe Trade park
Gateway
Crewe
CW1 6JT
UK

Tel: +44 (0)1270 580008
Reference: IN1294

Consignee:

LabTech Corporation, Inc.
No. 6 Anqing Street, Area B
Airport Industrial Zone
Shunyi District
Beijing, 101300
P.R. China

Contact Name: Elena Fan
Contact Tel Number: 010-80486450/1/2/3-8661

BOX: 1

Quantity	Product Code	Product Description	Serial Number: if available
☒ 1	20802/4	INTEGRA 2	20802/4-2023/1294

Packed by:

Date:

James Brown
27/1/23

Claims for damage or shortages must be made within 3 days



S035b_v6

Sercon Limited
Unit 3B Crewe Trade Park
Gateway, Crewe
CW1 6JT U.K.
Company Registration No.: 4290072

Phone: +44 (0) 1270 580008
Fax: +44 (0) 1270 252310
sales@sercongroup.com
www.sercongroup.com
VAT No.: 771 1100 74
19/01/2017



Packing List

Customer Reference: 20221000050

Shipper / Exporter:

Sercon Limited
Unit 3B Crewe Trade park
Gateway
Crewe
CW1 6JT
UK

Tel: +44 (0)1270 580008
Reference: IN1294

Consignee:

LabTech Corporation, Inc.
No. 6 Anqing Street, Area B
Airport Industrial Zone
Shunyi District
Beijing, 101300
P.R. China

Contact Name: Elena Fan
Contact Tel Number: 010-80486450/1/2/3-8661

BOX: 2

Quantity	Product Code	Product Description	Serial Number: if available
✓ 1		AGILENT PUMP + MAINS LEAD + UG15 OIL + FILTER	JP2220U036
✓ 1	SI0101	MAGNET	21/12/22
✓ 1	SI2391	MONITOR + MAINS LEAD	3QMP8H3
✓ 1	SI0155	PC + MAINS LEAD + USB (SOFTWARE/MANUALS/BACKUP)	8V2DQR3
✓ 1	10202-S	AUTOSAMPLER	42666
✓ 1		AUTOSAMPLER CONNECTIONS	
✓ 1		UNIT MAINS LEAD	
✓ 1	70201-2	ANALYSER SPARES KIT	
✓ 1	70501-3	SAMPLE PREP KIT	
✓ 1	70401-1	START UP CONSUMABLES KIT	
✓ 1	70601-3	SULPHUR OPTION	
✓ 1	70501-12	SIMULTANEOUS CNS CONSUMABLES KIT	

Packed by:

Date:

[Signature]

27/1/23

Claims for damage or shortages must be made within 3 days



S035b_v6

Sercon Limited
Unit 3B Crewe Trade Park
Gateway, Crewe
CW1 6JT U.K.
Company Registration No.: 4290072

Phone: +44 (0) 1270 580008
Fax: +44 (0) 1270 252310
sales@sercongroup.com
www.sercongroup.com
VAT No.: 771 1100 74
19/01/2017



Packing List

Customer Reference: 20221000050

Shipper / Exporter:

Sercon Limited
Unit 3B Crewe Trade park
Gateway
Crewe
CW1 6JT
UK

Tel: +44 (0)1270 580008
Reference: IN1294

Consignee:

LabTech Corporation, Inc.
No. 6 Anqing Street, Area B
Airport Industrial Zone
Shunyi District
Beijing, 101300
P.R. China

Contact Name: Elena Fan
Contact Tel Number: 010-80486450/1/2/3-8661

BOX: Dimensions and weights

	<u>L</u>		<u>W</u>		<u>H</u>		
Box 1	100	x	86	x	110	cm	gw 154 kg
Box 2	100	x	86	x	110	cm	gw 80 kg



S035b_v6

Sercon Limited
Unit 3B Crewe Trade Park
Gateway, Crewe
CW1 6JT U.K.
Company Registration No.: 4290072

Phone: +44 (0) 1270 580008
Fax: +44 (0) 1270 252310
sales@sercongroup.com
www.sercongroup.com
VAT No.: 771 1100 74
19/01/2017



COLLECTION NOTE

Job Number: **IN1294**
Customer: Weifang Nongfeng Testing Co., Ltd C/O: LabTech
Contact: **Elena Fan**

Collection on
Date:

2/2/23

No.
of boxes:

2

By:

Company:

UK FREIGHT & LOGISTICS

Vehicle

Reg No.

DV21VJA

Signature:

[Signature]

Print:

S. SINGH



S062_v3

Sercon Ltd
Unit 3B Crewe Trade Park
Gateway, Crewe, Cheshire
CW1 6JT U.K.
Registered in England and Wales: 4290072

Phone: +44 (0) 1270 580008
Fax: +44 (0) 1270 252310
sales@sercongroup.com
www.sercongroup.com
VAT No.: 771 1100 74
19/01/2017